




# Inspired by temperature

Betriebsanleitung · Operation manual · Manual de instrucciones · Manuel d'utilisation · Manuale de d'uso · 사용 설명서 · Manual de instruções · Инструкция по эксплуатации · Kullanım talimatı · 操作说明书 · Betriebsanweisung · Manual de instrucciones · 사용 설명서 · Инструкция по эксплуатации · Kullanım talimatı · 操作说明书 · Betriebsanweisung · Manual de instrucciones · 사용 설명서 · Инструкция по эксплуатации · Kullanım talimatı · 操作说明书

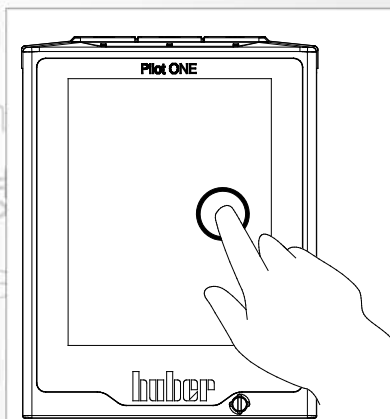

**CC<sup>®</sup>**  
Immersion  
Circulator

**This documentation does not contain a device-specific technical appendix.**


You can request the full operating instructions from [info@huber-online.com](mailto:info@huber-online.com). Please give the model designation and serial number of your temperature control unit in your e-mail.

**huber**



# huber

## OPERATION MANUAL

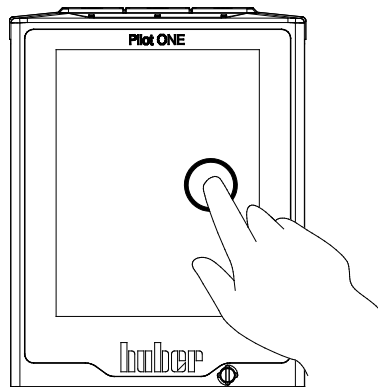

**CC<sup>®</sup>**  
Immersion  
Circulator



# CC<sup>®</sup>-Immersion circulator

Pilot ONE<sup>®</sup>

This operation manual is a translation of the original operation manual.


## VALID FOR:

CC<sup>®</sup>-E  
CC<sup>®</sup>-1xx  
CC<sup>®</sup>-130 Visco x  
CC<sup>®</sup>-2xx  
CC<sup>®</sup>-3xx  
CC<sup>®</sup>-2xxx

Abbreviations used in model names:


A = polycarbonate bath, B = stainless steel bath, BX = telescope arm, C = cooling coil,  
xd = greater immersion depth

„Pilot ONE“

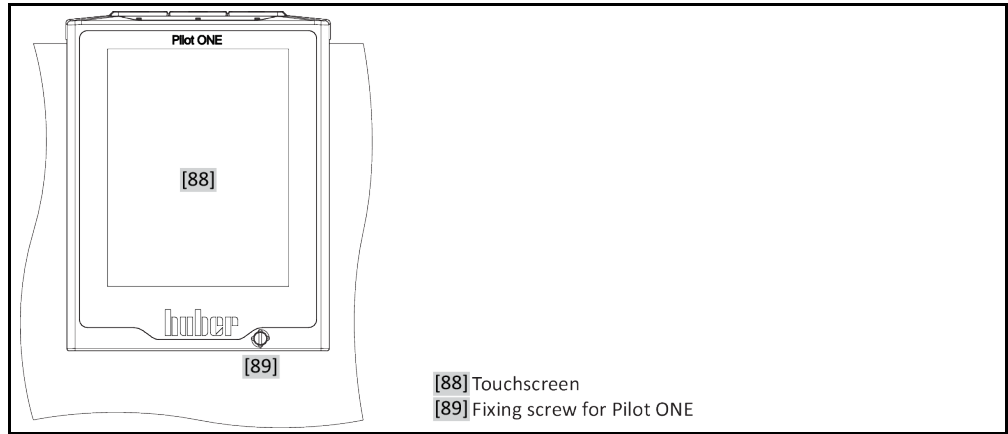

- [88] Touchscreen
- [89] Fixing screw for Pilot ONE

Layout of the “Home” screen

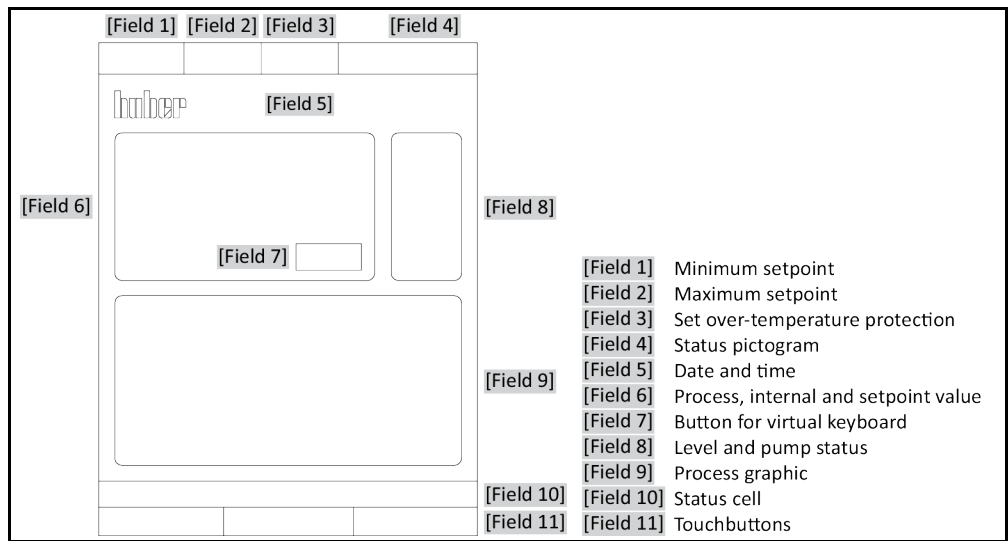

- [Field 1] Minimum setpoint
- [Field 2] Maximum setpoint
- [Field 3] Set over-temperature protection
- [Field 4] Status pictogram
- [Field 5] Date and time
- [Field 6] Process, internal and setpoint value
- [Field 7] Button for virtual keyboard
- [Field 8] Level and pump status
- [Field 9] Process graphic
- [Field 10] Status cell
- [Field 11] Touchbuttons

# Table of contents

V2.4.0en/08.12.20//17.12

<b>1</b>	<b>Introduction</b>	<b>12</b>
<b>1.1</b>	<b>Details on the declaration of conformity.....</b>	<b>12</b>
<b>1.2</b>	<b>Safety .....</b>	<b>12</b>
1.2.1	Symbols used for Safety Instructions.....	12
1.2.2	Representation of safety identifiers .....	13
1.2.3	Proper operation.....	13
1.2.4	Reasonably foreseeable misuse.....	14
<b>1.3</b>	<b>Responsible bodies and operators – Obligations and requirements.....</b>	<b>14</b>
1.3.1	Obligations of the responsible body .....	14
1.3.1.1	Proper disposal of resources and consumables .....	15
1.3.2	Requirements for operators .....	15
1.3.3	Obligations of the operators.....	15
<b>1.4</b>	<b>General information.....</b>	<b>15</b>
1.4.1	Description of workstation .....	15
1.4.2	Safety devices to DIN 12876 .....	16
1.4.3	Further protective devices.....	17
1.4.3.1	Power interruption.....	17
1.4.3.2	Alarm functions .....	17
1.4.3.3	Warning messages.....	17
<b>2</b>	<b>Commissioning</b>	<b>18</b>
<b>2.1</b>	<b>In-plant transport.....</b>	<b>18</b>
2.1.1	Lifting and transporting the temperature control unit.....	18
2.1.1.1	Temperature control unit with lifting eyes .....	18
2.1.1.2	Temperature control unit without lifting eyes.....	19
2.1.2	Mounting/removing leveling feet.....	19
2.1.3	Positioning the temperature control unit.....	20
2.1.3.1	Temperature control unit with casters .....	20
2.1.3.2	Temperature control unit without casters.....	20
<b>2.2</b>	<b>Unpacking.....</b>	<b>20</b>
<b>2.3</b>	<b>Ambient conditions.....</b>	<b>20</b>
2.3.1	EMC-specific notes.....	22
<b>2.4</b>	<b>Installation conditions.....</b>	<b>22</b>
<b>2.5</b>	<b>Recommended temperature control and cooling water hoses .....</b>	<b>23</b>
<b>2.6</b>	<b>Wrench sizes and torques.....</b>	<b>23</b>
<b>2.7</b>	<b>Temperature control units with cooler .....</b>	<b>24</b>
<b>2.8</b>	<b>Temperature control unit with batch.....</b>	<b>25</b>
2.8.1	Operation as bath thermostat .....	25
<b>2.9</b>	<b>Preparations for operation .....</b>	<b>25</b>
2.9.1	Unscrewing/activating the leveling feet (if any) .....	25
2.9.2	Externally closed and externally open applications.....	25
2.9.3	Opening/closing valves .....	26
2.9.4	Installing collecting container .....	26
2.9.5	Connecting the functional earth .....	26
<b>2.10</b>	<b>Connecting externally closed application.....</b>	<b>26</b>
2.10.1	Connecting an externally closed application .....	26
<b>2.11</b>	<b>Connecting to the power supply .....</b>	<b>27</b>

2.11.1	Connection using socket with protective earth (PE).....	27
2.11.2	Connection via hard wiring.....	27
2.11.3	Converting the power supply connection.....	28
<b>3</b>	<b>Function description</b>	<b>29</b>
<b>3.1</b>	<b>Function description of the temperature control unit</b> .....	<b>29</b>
3.1.1	General functions.....	29
3.1.2	Other functions.....	29
<b>3.2</b>	<b>Information on the thermal fluids</b> .....	<b>29</b>
<b>3.3</b>	<b>To be noted when planning the test</b> .....	<b>30</b>
<b>3.4</b>	<b>“Pilot ONE®” controller</b> .....	<b>31</b>
3.4.1	Functional overview of “Pilot ONE®” .....	31
<b>3.5</b>	<b>Clock/event function</b> .....	<b>33</b>
3.5.1	Rechargeable accumulator .....	33
3.5.2	Programmable event function .....	33
3.5.2.1	Event function “Alarm clock event” .....	33
3.5.2.2	Event function “Program event” .....	33
<b>3.6</b>	<b>Operation via the touch screen</b> .....	<b>34</b>
<b>3.7</b>	<b>Display instruments</b> .....	<b>34</b>
3.7.1	The touchscreen [88] .....	34
<b>3.8</b>	<b>Control instruments</b> .....	<b>34</b>
3.8.1	The touchbuttons.....	35
3.8.2	The categories.....	35
3.8.3	The sub-categories.....	35
3.8.4	The dialogs .....	35
<b>3.9</b>	<b>Function examples</b> .....	<b>35</b>
3.9.1	Display of software version.....	35
3.9.2	Start & Stop.....	36
3.9.3	Copying the settings to a data carrier.....	36
3.9.3.1	Saving to a USB flash drive .....	37
3.9.3.2	Loading from a USB flash drive.....	37
3.9.4	Restore factory settings .....	37
3.9.4.1	Restore to factory settings without overtemperature protection .....	38
3.9.4.2	Restore to factory settings including overtemperature protection .....	39
<b>4</b>	<b>Setup mode</b>	<b>40</b>
<b>4.1</b>	<b>Setup mode</b> .....	<b>40</b>
4.1.1	Turning on the temperature control unit .....	40
4.1.2	Turning off the temperature control unit.....	40
4.1.3	Setting the overtemperature protection .....	41
4.1.3.1	General information on the overtemperature protection .....	41
4.1.3.2	Setting “OT limit: heating” .....	41
4.1.3.3	Setting “Process Safety” .....	42
4.1.3.4	Monitoring via “Display OT values” .....	42
4.1.4	Testing overtemperature protection for functionality .....	42
4.1.5	Adjusting the Delta T limiter .....	43
4.1.5.1	Changing the Delta T limiter.....	43
<b>4.2</b>	<b>The temperature control circuit</b> .....	<b>43</b>
4.2.1	Select temperature control: Internal or process .....	44
4.2.2	Temperature control to internal temperature .....	44
4.2.3	Temperature control to process temperature.....	44
4.2.4	Delta T limiter .....	45
4.2.5	Monitoring the Pt100 temperature sensors .....	45



4.2.6	Optimum control parameters for optimum temperature control.....	45
4.2.7	Sub-category: "Select auto/expert mode" .....	46
4.2.8	Sub-category: "Configuration auto" .....	46
4.2.8.1	Sub-category: "Find parameters" .....	46
4.2.8.2	Sub-category: "Control Dynamics" .....	48
4.2.8.3	Sub-category: "Fluid Properties" .....	49
4.2.8.4	Sub-category: "Display parameters" .....	50
4.2.9	Sub-category: "Configuration manual" .....	50
4.2.9.1	Sub-category: "Change parameters" .....	50
4.2.9.2	Sub-category: "Display parameters" .....	51
4.2.9.3	Sub-category: "Control structure" .....	52
4.2.10	Sub-category: "Reset parameters" .....	52
4.2.11	Sub-category: "Display parameters" .....	52
4.2.12	Freeze protection for temperature control unit .....	52
4.2.13	Setting the setpoint thresholds .....	53
4.2.14	Setting the setpoint .....	54
<b>4.3</b>	<b>Filling, venting, degassing and draining.....</b>	<b>54</b>
4.3.1	Filling, venting, degassing and draining the bath thermostat.....	54
4.3.1.1	Filling and venting the bath thermostat and the externally closed application .....	55
4.3.1.2	Degassing of bath thermostat .....	55
4.3.1.3	Draining the bath thermostat.....	56
<b>5</b>	<b>Normal operation .....</b>	<b>58</b>
<b>5.1</b>	<b>Automatic operation .....</b>	<b>58</b>
5.1.1	Temperature control.....	58
5.1.1.1	Starting the temperature control process.....	58
5.1.1.2	Ending the temperature control process .....	58
5.1.2	Temperature control via a created temperature control program .....	59
5.1.2.1	Starting the temperature control program .....	59
5.1.2.2	Ending/cancelling the temperature control program .....	59
<b>6</b>	<b>Interfaces and software update .....</b>	<b>60</b>
<b>6.1</b>	<b>Interfaces at the "Pilot ONE®" controller .....</b>	<b>60</b>
6.1.1	10/100 Mbps Ethernet for RJ45 network sockets .....	60
6.1.2	USB-2.0 interface .....	61
6.1.2.1	USB-2.0 interface, host.....	61
6.1.2.2	USB-2.0 interface, device .....	61
<b>6.2</b>	<b>Interfaces on the temperature control unit.....</b>	<b>61</b>
6.2.1	Interfaces at the back .....	61
6.2.1.1	Service interface.....	61
6.2.1.2	Female RS232 serial.....	61
6.2.1.3	Connection jack for Pt100 process controller sensor .....	62
<b>6.3</b>	<b>Firmware update.....</b>	<b>62</b>
<b>7</b>	<b>Service/maintenance .....</b>	<b>63</b>
<b>7.1</b>	<b>Messages from the temperature control unit .....</b>	<b>63</b>
<b>7.2</b>	<b>Electrical fuse.....</b>	<b>63</b>
<b>7.3</b>	<b>Replacing the "Pilot ONE®" .....</b>	<b>63</b>
<b>7.4</b>	<b>Maintenance.....</b>	<b>64</b>
7.4.1	Function check and visual inspection .....	64
7.4.2	Replacing temperature control or coolant hoses .....	65
7.4.2.1	Replacing temperature control hoses .....	65

7.4.2.2	Replacing coolant hoses .....	65
<b>7.5</b>	<b>Thermal fluid inspection, replacement and circuit cleaning .....</b>	<b>65</b>
7.5.1	Thermal fluid inspection .....	65
7.5.2	Rinsing the thermal fluid circuit .....	66
<b>7.6</b>	<b>Cleaning the surfaces .....</b>	<b>67</b>
<b>7.7</b>	<b>Plug contacts.....</b>	<b>68</b>
<b>7.8</b>	<b>Decontamination/repairs .....</b>	<b>68</b>
<b>8</b>	<b>Shutting down .....</b>	<b>69</b>
<b>8.1</b>	<b>Safety instructions and basic principles .....</b>	<b>69</b>
<b>8.2</b>	<b>Switch-off .....</b>	<b>69</b>
<b>8.3</b>	<b>Draining the temperature control unit.....</b>	<b>70</b>
<b>8.4</b>	<b>Draining the cooling water .....</b>	<b>70</b>
8.4.1	Draining process .....	70
<b>8.5</b>	<b>Deinstalling the collecting container .....</b>	<b>70</b>
<b>8.6</b>	<b>Uninstalling an external application .....</b>	<b>70</b>
<b>8.7</b>	<b>Packing .....</b>	<b>70</b>
<b>8.8</b>	<b>Shipping.....</b>	<b>71</b>
<b>8.9</b>	<b>Disposal .....</b>	<b>71</b>
<b>8.10</b>	<b>Contact data .....</b>	<b>72</b>
8.10.1	Telephone number: Customer Support .....	72
8.10.2	Telephone number: Sales .....	72
8.10.3	Email address: Customer Support.....	72
<b>8.11</b>	<b>Certificate of Compliance .....</b>	<b>72</b>
<b>9</b>	<b>Annex .....</b>	<b>73</b>

## Foreword

Dear Customer,

Thank you for choosing a temperature control unit from Peter Huber Kältemaschinenbau AG. You have made a good choice. Thank you for your trust.

Please read the operation manual carefully before putting the unit into operation. Strictly follow all notes and safety instructions.

Follow the operation manual with regard to transport, start-up, operation, maintenance, repair, storage and disposal of the temperature control unit.

We fully warrant the temperature control unit for the specified intended operation.


The models listed on page 5 are referred to in this operation manual as temperature control units and Peter Huber Kältemaschinenbau AG as Huber company or Huber.

Liability for errors and misprints excluded.

The following trademarks and the Huber logo are registered trademarks of Peter Huber Kältemaschinenbau AG in Germany and/or other countries worldwide: BFT®, CC®, Chili®, Com.G@te®, Compatible Control®, CoolNet®, DC®, E-grade®, Grande Fleur®, KISS®, Minichiller®, Ministat®, MP®, MPC®, Peter Huber Minichiller®, Petite Fleur®, Pilot ONE®, RotaCool®, Rotostat®, SpyControl®, SpyLight®, Tango®, TC®, UC®, Unical®, Unichiller®, Unimotive®, Unipump®, Unistat®, Unistat Tango®, Variostat®. The following trademarks are registered in Germany to DWS Synthesetechnik: DW-Therm®, DW-Therm HT®. The following trademark is a registered trademark of BASF SE: Glysantin®.

# 1 Introduction

## 1.1 Details on the declaration of conformity





 The equipment complies with the basic health and safety requirements of the European Directives listed below:

- Machinery Directive
- Low Voltage Directive
- EMC Directive


## 1.2 Safety

### 1.2.1 Symbols used for Safety Instructions

Safety instructions are marked by the below combinations of pictograms and signal words. The signal word describes the classification of the residual risk when disregarding the operation manual.

 <b>DANGER</b>	Denotes an immediate hazardous situation that will result in death or serious injuries.
 <b>WARNING</b>	Denotes a general hazardous situation that may result in death or serious injuries.
 <b>CAUTION</b>	Denotes a hazardous situation that can result in injury.
<b>NOTE</b>	Denotes a situation that can result in property material damage.
<b>INFORMATION</b>	Denotes important notes and usable hints.
	Notes in conjunction with Ex px cabinets.

Safety information and procedure

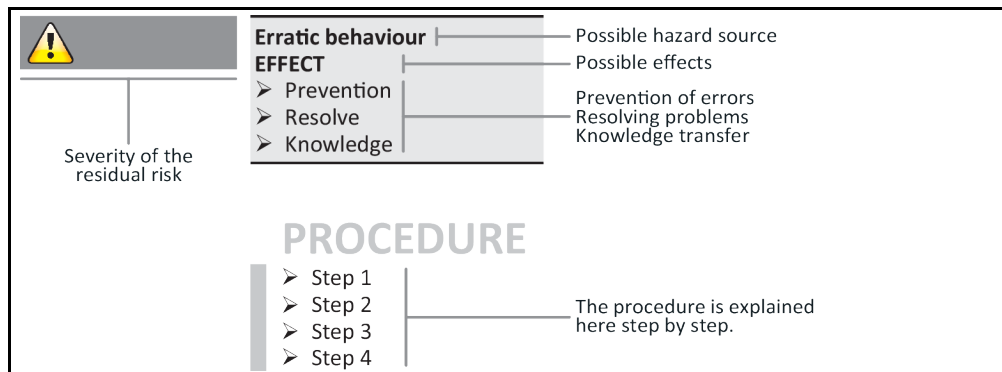

The safety information in this operation manual is designed to protect the responsible body, the operator and the equipment from damage. Safety instructions appear at the beginning of each chapter and before instructions. First inform yourself about any residual risks due to misuse before you start an operation.

### 1.2.2 Representation of safety identifiers

The following pictograms are used as safety identifiers. The table gives an overview of the safety identifier used here.

Identifier	Description
<b>Mandatory sign</b>	
	- Observe the instructions
<b>Warning sign</b>	
	- General warning sign - observe the instructions
	- Warning of electrical voltage
	- Warning of hot surface
	- Warning of flammable substances

### 1.2.3 Proper operation



**Operating the temperature control unit in a potentially explosive area**  
**DEATH THROUGH EXPLOSION**  
 ➤ Do NOT install or start up the temperature control unit within an ATEX zone.



**Improper use**  
**SERIOUS INJURY AND PROPERTY DAMAGE**  
 ➤ Store the operation manual where it is easy to access in close proximity to the temperature control unit.  
 ➤ Only adequately qualified operators may work with the temperature control unit.  
 ➤ Operators must be trained before handling the temperature control unit.  
 ➤ Check that the operators have read and understood the operation manual.  
 ➤ Define precise responsibilities of the operators.  
 ➤ Personal protective equipment must be provided to the operators.  
 ➤ Be sure to follow the responsible body's safety rules to protect life and limb and to limit damages!



**Modifications to the temperature control unit by third-parties**  
**DAMAGE TO THE TEMPERATURE CONTROL UNIT**  
 ➤ Do not allow third parties to make technical modifications to the temperature control unit.  
 ➤ The EU declaration of conformity becomes invalid if any modifications are made to the temperature control unit without the approval of Huber.  
 ➤ Only specialists trained by Huber may carry out modifications, repairs or maintenance work.  
 ➤ **The following must be observed without fail:**  
 ➤ Only use the temperature control unit in a fault-free condition!  
 ➤ Have the start-up and repairs carried out by specialists only!  
 ➤ Do not ignore, bypass, dismantle or disconnect any safety devices!

The temperature control unit must not be used for any purposes other than temperature control in accordance with the operation manual.

The temperature control unit is manufactured for industrial use. The temperature control unit maintains the temperature of certain applications, including glass or metal reactors or other expedient items in laboratories and industry. Flow-through coolers and calibration baths must be used only in combination with Huber temperature control units. Thermal fluids suitable for the overall system are used. The chilling and heating capacity is provisioned at the pump connections or - where present - in the tempering bath. The technical specification of the temperature control unit is given in the data sheet (from page 73 in section »Annex«). The temperature control unit must be installed, configured and operated according to the handling instructions in this operation manual. Failure to comply with the operation manual is deemed improper use. The temperature control unit conforms to state-of-the-art technology and the recognized safety regulations. Safety devices are built into your temperature control unit.

#### 1.2.4 Reasonably foreseeable misuse



Without an Ex px cabinet, the temperature control unit / accessory is **NOT** protected against explosion and must **NOT** be installed or put into operation within an ATEX Zone. When operating the temperature control unit /accessory in conjunction with an Ex px cabinet, the information in the annex (Section ATEX operation) must be observed and followed. This annex is only provided for temperature control units /accessories delivered with an Ex px cabinet. If this annex is missing, please immediately contact the Customer Support of Huber (the telephone number is provided on page 72 in section »Contact data«).

Use with medical devices (e.g. in Vitro diagnostic procedure) or for direct foodstuff temperature control is **NOT** permissible.

The temperature control unit must **NOT** be used for any purposes other than temperature control in accordance with the operation manual.

The manufacturer accepts **NO** liability for damage caused by **technical modifications** to the temperature control unit, **improper handling** or use of the temperature control unit if the operation manual is **not observed**.

### 1.3 Responsible bodies and operators – Obligations and requirements

#### 1.3.1 Obligations of the responsible body

The operation manual is to be stored where it is easy to access in close proximity to the temperature control unit. Only adequately qualified operators (e.g. chemists, CTA, physicists etc.) are permitted to work with the temperature control unit. Operators must be trained before handling the temperature control unit. Check that the operators have read and understood the operation manual. Define precise responsibilities of the operators. Personal protective equipment must be provided to the operators.

- The responsible body must install a condensation water / thermal fluid drip tray below the temperature control unit.
- The use of a drain tray may be prescribed by national law for the installation area of the temperature control unit (including accessories). The responsible body must check and apply the applicable national regulations.
- Our temperature control unit complies with all applicable safety standards.
- Your system, which uses our temperature control unit, must be as safe.
- The responsible body must design the system to ensure it is safe.
- Huber is not responsible for the safety of your system. The responsible body is responsible for the safety of the system.
- Although the temperature control unit provided by Huber meets all the applicable safety standards, integration into a system may give rise to hazards that are characteristic of the other sys-

tem’s design and beyond the control of Huber.

- It is the responsibility of the system integrator to ensure that the overall system, into which this temperature control unit is integrated, is safe.
- The >Mains isolator< [36] (if present) may be provided with a facility to lock the main isolator in the off position to facilitate safe system installation and maintenance of the temperature control unit. It is the responsibility of the responsible body to develop any lock-out/tag-out procedure in accordance with local regulations (e.g. CFR 1910.147 for the US).

**1.3.1.1 Proper disposal of resources and consumables**

Do comply with all national disposal regulations applicable for you. Contact your local waste management company for any questions concerning disposal.

Material / Aids	Disposal / Cleaning
Packing material	Keep the packaging material for future use (e.g. transport).
Thermal fluid	Please refer to the safety data sheet of the thermal fluid used for information on its proper disposal. Use the original thermal fluid container when disposing it.
Filling accessories, e.g. beaker	Clean the filling accessories for reuse. Make sure that the materials and cleaning agents used are properly disposed of.
Aids such as towels, cleaning cloths	Tools used to take up spilled thermal fluid must be disposed of in the same fashion as the thermal fluid itself. Tools used for cleaning must be disposed of depending on the cleaning agent used.
Cleaning agents such as stainless steel cleaning agents, sensitive-fabrics detergents	Please refer to the safety data sheet of the cleaning agent used for information on its proper disposal. Use the original containers when disposing of large quantities of cleaning agents.
Consumables such as air filter mats, temperature control hoses	Please refer to the safety data sheet of the consumables used for information on their proper disposal.

**1.3.2 Requirements for operators**

Work on the temperature control unit is reserved for appropriately qualified specialists, who have been assigned and trained by the responsible body to do so. Operators must be at least 18 years old. Under 18-year olds may operate the temperature control unit only under the supervision of a qualified specialist. The operator is responsible vis-a-vis third-parties in the work area.

**1.3.3 Obligations of the operators**

Carefully read the operation manual before operating the temperature control unit. Please observe the safety instructions. When operating the temperature control unit, wear appropriate personal protective equipment (e.g. safety goggles, protective gloves, non-slip shoes).

**1.4 General information**

**1.4.1 Description of workstation**

The workstation is located at the control panel in front of the temperature control unit. The workstation is determined by the customer's connected peripheries. Accordingly, it must be designed safe by the responsible body. The workstation design also depends on the applicable requirements of the German occupational health and safety regulations [BetrSichV] and the risk analysis for the workstation.


### 1.4.2 Safety devices to DIN 12876

The rating of your temperature control unit is stated on the data sheet in the appendix.

Rating of laboratory thermostats and laboratory baths	Classification	Temperature control medium	Technical requirements	Identification <sup>d)</sup>
	I	Non-combustible <sup>a)</sup>	Overheat protection <sup>c)</sup>	NFL
	II	Combustible <sup>b)</sup>	Adjustable overheat protection	FL
	III	Combustible <sup>b)</sup>	Adjustable overtemperature protection and additional low-level protection	FL

<sup>a)</sup> Usually water; other fluids only if non-combustible even within the temperature range of an individual fault.  
<sup>b)</sup> The temperature control media must have a fire point of  $\geq 65$  °C.  
<sup>c)</sup> The overheat protection can, for instance, can be realized using a suitable fill level sensor or a suitable temperature limiter.  
<sup>d)</sup> Optional at the choice of the manufacturer.

Overview of the temperature thresholds

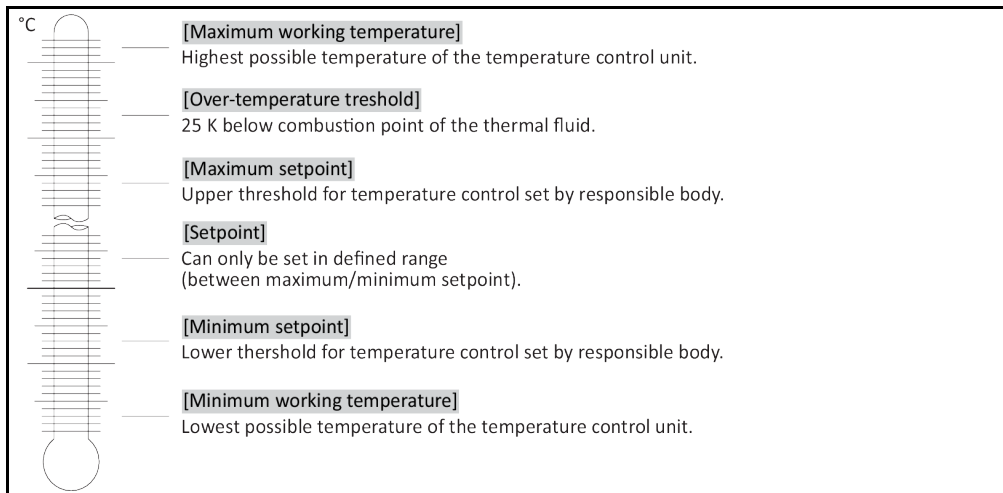

#### Low level protection

Probably the most frequently found type today is the level monitor with a mechanical float switch. In the bath vessel, a floating body, which is guided in a device, floats on the surface of the thermal fluid. Depending on the level of the thermal fluid, the float device signals the electronics a **state of good** (in case of sufficient filling) or a **state of bad** (in case of insufficient filling). You should check the functionality of the float device from time to time. To do so and when in stand-by mode, press the float body into the bath with a tool (e.g. a screwdriver). The electronics must report an alarm.

#### Overtemperature protection (for temperature control units with heating)

Mechanical tools are no longer required for setting the trip values for the overtemperature protection. A software tool is used instead. The threshold value for the overtemperature protection can be set only if a code randomly generated by "Pilot ONE" beforehand is entered correctly. As with the mechanical tool, accidental settings are thus prevented. The type of overtemperature protection function and of low level protection depends on the temperature control unit.



### 1.4.3 Further protective devices

#### INFORMATION

Emergency strategy – isolate the power supply!

To determine the type of switch or switch combination your temperature control unit is equipped with, please refer to the connection sketch starting on page 73 in the section »Annex«.

**Temperature control units with >Mains isolator< [36] (red/yellow or gray):** Turn the >Mains isolator< [36] to the "0" position!

**Temperature control units with >Mains isolator< [36] (red/yellow) and additional >Appliance switch< [37] (gray):** Turn the >Mains isolator< [36] to the "0" position! Then turn the >Appliance switch< [37] to the "0" position!

**Temperature control units with >Mains isolator< [36] (gray) and >Emergency stop switch< [70] (red/yellow):** Press the >Emergency stop switch< [70]. Then turn the >Mains isolator< [36] to the "0" position!

**Temperature control units with >Mains switch< [37]:** Power supply via socket: Disconnect the temperature control unit from the power supply. Then turn the >Mains isolator< [37] to the "0" position! Power supply via hard wiring: Disconnect the power supply by means of the building's circuit breaker. Then turn the >Mains isolator< [37] to the "0" position!

**Temperature control units without a switch or inside a protective housing:** Connection via socket: Disconnect the temperature control unit from the power supply. Connection via hard wiring: Disconnect the power supply by means of the building's circuit breaker.

#### 1.4.3.1 Power interruption

Following a power outage (or when switching on the temperature control unit), this function can be used to determine how the temperature control unit is supposed to respond. This response can be determined via "Pilot ONE".

**OFF/Standby** (Default setting)

After turning the temperature control unit on, thermoregulation is started only after manual input.

**ON / Temperature control active**

After turning on the temperature control unit, thermoregulation is always started. An INFO appears for a few seconds. This makes it possible to suppress the automatic start.

**Power Failure Auto**

If temperature control is active during a power outage, the process will automatically continue after the power outage.

#### 1.4.3.2 Alarm functions

An alarm is a system state that signals unfavorable process conditions. The temperature control unit can be programmed so that the plant operator is warned when defined limit values are exceeded.

The response of the temperature control unit to an alarm can be determined. Possible responses are: Switch off temperature control or control temperature to a safe setpoint (2nd setpoint).

#### 1.4.3.3 Warning messages

Warning messages contain a message about the irregularity of the temperature control unit. These messages have no further consequences. The plant operator evaluates the relevance of the message and takes action where necessary.

## 2 Commissioning

### 2.1 In-plant transport

**WARNING**

Temperature control unit is not transported / moved according to the specifications in this operation manual

**DEATH OR SERIOUS INJURY DUE TO CRUSHING**

- Always transport / move the temperature control unit according to the specifications in this operation manual.
- Wear personal protective equipment during transport.
- Always work with the specified number of persons when moving the temperature control unit on casters (if any).
- If the temperature control unit is equipped with casters and parking brakes: 2 parking brakes are always freely accessible when moving the temperature control unit. Activate the **2 parking brakes** in an emergency! If only **one** parking brake is activated on the casters in an emergency: The temperature control unit is not stopped but rotates around the axis of the caster with the activated parking brake!

**NOTE**

Temperature control unit transported in a horizontal position

**DAMAGE TO THE COMPRESSOR**

- Only transport the temperature control unit in an upright position.

**NOTE**

Filled temperature control unit is transported

**MATERIAL DAMAGE DUE TO OVERFLOWING THERMAL FLUID**

- Only transport an emptied temperature control unit.
- If available, use the lugs on the top side of the temperature control unit for transportation.
- Use an industrial truck for transport.
- The casters (if present) on the temperature control unit are not suitable for transport. The casters are symmetrically loaded with 25% of the total mass of the temperature control unit.
- Remove the packing material (e.g. the palette) only at the place of installation.
- Protect the temperature control unit from transport damage.
- Do not transport the temperature control unit alone and without aids.
- Check the load bearing capacity of the transportation route and the place of installation.
- The parking brakes must be activated at the casters (if any) and/or the leveling feet (if any) must be unscrewed/activated before the temperature control unit is put into operation (see page 25, section »Unscrewing/activating the leveling feet (if any)«).

#### 2.1.1 Lifting and transporting the temperature control unit

##### 2.1.1.1 Temperature control unit with lifting eyes


**NOTE**

The temperature control unit is raised at the lifting eyes without load handling attachments

**DAMAGE TO THE TEMPERATURE CONTROL UNIT**

- Always use load handling attachments when lifting and transporting the temperature control unit.
- The lifting eyes are only designed for a load **without** inclination (0°).
- The load handling attachment used must be adequately dimensioned. Take the dimensions and weight of the temperature control unit into account.


Example: lifting eyes (round, angular, and recessed (left to right))

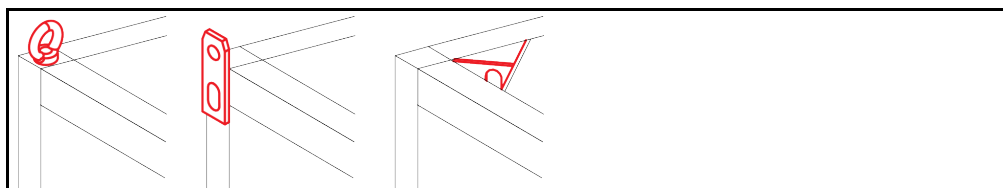

- Do not lift and transport the temperature control unit at the lifting eyes alone and without aids.
- Lift and transport the temperature control unit at the lifting eyes only with a crane or an industrial truck.
- The crane or industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. For the weight of the temperature control unit see the data sheet (from page 73 in section »Annex«).
- If the leveling feet were removed for shipment: Do not lower the temperature control unit until all leveling feet have been mounted (see page 19 section »Mounting/removing leveling feet«).

**2.1.1.2 Temperature control unit without lifting eyes**

Example: Supporting points for forklift arms for free-standing models from a certain overall size. For the exact position please refer to the wiring diagram in the annex.

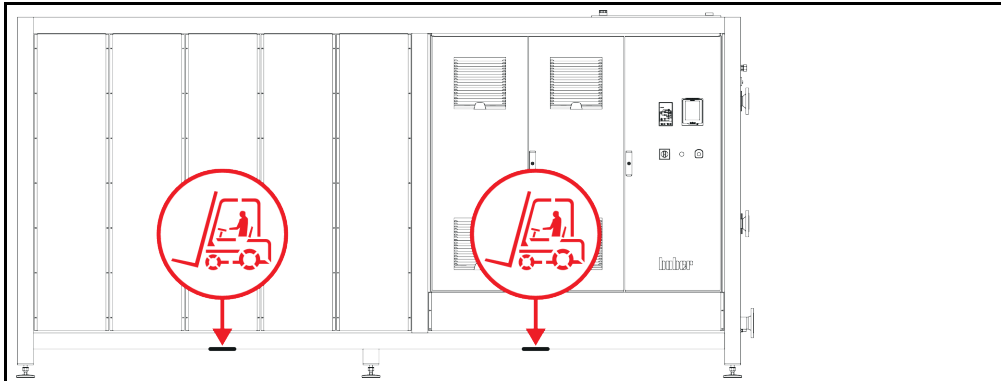

- Do not lift and transport the temperature control unit alone and without aids.
- Lift and transport the temperature control unit only with an industrial truck.
- The industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. For the weight of the temperature control unit see the data sheet (from page 73 in section »Annex«).
- If the leveling feet were removed for shipment: Do not lower the temperature control unit until all leveling feet have been mounted (see page 19 section »Mounting/removing leveling feet«).

**2.1.2 Mounting/removing leveling feet**

Only valid if the leveling feet have been removed for shipping.




**WARNING**

**The temperature control unit is not secured against slipping and/or lowering  
DEATH OR SERIOUS INJURY DUE TO CRUSHING**

- Secure the temperature control unit against slipping and/or lowering before the leveling feet are mounted.
- Do not stand or lie under the temperature control unit for mounting.

Example: mounting the leveling feet

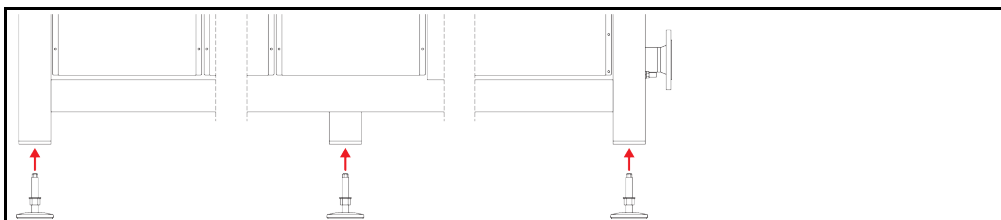

**INFORMATION**

The leveling feet were removed for shipping the temperature control unit. Before placing / positioning the temperature control unit all leveling feet must be mounted. If the temperature control unit is re-shipped: Remove all leveling feet before packaging.

- The leveling feet can only be mounted while the temperature control unit is lifted.
- Secure the temperature control unit against slipping and/or lowering.
- Do not stand or lie under the temperature control unit while mounting the leveling feet.
- Do not lower the temperature control unit until all leveling feet have been mounted.

## 2.1.3 Positioning the temperature control unit

### 2.1.3.1 Temperature control unit with casters

- Do **not** use the casters for transportation to the place of installation. Observe page 18, section »**Lifting and transporting the temperature control unit**« for the transport to the place of installation.
- Use the rollers only for positioning at the place of installation.
- Only ever move the temperature control unit on casters if the surface is level, without a gradient, non-slip and stable.
- Do not move the temperature control unit alone.
- **At least 2 persons** are required to move the temperature control unit on casters. **At least 5 persons** are required to move the temperature control unit on casters if the total weight of the temperature control unit is **over 1.5 tons**.
- The parking brakes must be activated at the casters and/or the leveling feet (if any) must be unscrewed/activated before the temperature control unit is put into operation (see page 25, section »**Unscrewing/activating the leveling feet (if any)**«).

### 2.1.3.2 Temperature control unit without casters

- An industrial truck must be used for positioning the temperature control unit.
- Do not move the temperature control unit alone.
- **At least 2 persons** are required to move the temperature control unit.
- The industrial truck must have a lifting force equal to or greater than the weight of the temperature control unit. See the data sheet (from page 73 in section »**Annex**«) for the weight of the temperature control unit.
- The leveling feet (if any) must be unscrewed/activated before the temperature control unit is put into operation (see page 25, section »**Unscrewing/activating the leveling feet (if any)**«).

## 2.2 Unpacking

### WARNING

#### Starting up a damaged temperature control unit

#### MORTAL DANGER FROM ELECTRIC SHOCK

- Do not operate a damaged temperature control unit.
- Please contact the Customer Support. The telephone number can be found on page 72, section »**Contact data**«.

## PROCEDURE

- Check for damage to the packaging. Damage can indicate property damage to the temperature control unit.
- Check for any transport damage when unpacking the temperature control unit.
- Always contact your forwarding agent regarding the settlement of claims.
- Follow the instructions for the disposal of packaging material on page 15 section »**Proper disposal of resources and consumables**«.

## 2.3 Ambient conditions

### CAUTION

#### Unsuitable ambient conditions/unsuitable installation

#### SERIOUS INJURY DUE TO CRUSHING

- Comply with the requirements under sections »**Ambient conditions**« and »**Installation conditions**«.

### INFORMATION

Make sure there is adequate fresh air available at the site for the circulation pump and the compressors. The warm exhaust air must be able to escape upwards unhindered.

#### Free-standing model

For the connection data, see the data sheet (from page 73 in Section »**Annex**«).

Use of the temperature control unit is permitted only under normal ambient conditions in accordance with DIN EN 61010-1:2011:

- Use only indoors. The illuminance must be at least 300 lx.
- Installation altitude up to 2000 meters above sea level .
- Maintain wall and ceiling clearance for adequate air exchange (dissipation of waste heat, supply of fresh air for the temperature control unit and work area). Ensure adequate floor clearance for air-cooled temperature control units. Do not operate this temperature control unit from within the box or with an inadequately dimensioned bath. This inhibits the exchange of air.
- Ambient temperature values are provided on the technical data sheet; compliance with the ambient conditions is mandatory, to ensure trouble-free operation.
- Relative humidity up to 32 °C max. 80% and decreasing linearly to 50% up to 40 °C.
- Short distance to supply connections.
- The temperature control unit must not be installed so as to hinder or prevent access to the isolator (to the power supply).
- Magnitude of the power supply fluctuations: see data sheet from page 73 in section »Annex«.
- Transient surges, as would normally occur in the power supply system
- Installation Class 3
- Applicable pollution degree: 2.
- Surge category II.

Wall clearances

Distance in cm	
<b>Side</b>	
[A2] Top	can be located under a bench
[B] Left	min. 10
[C] Right	min. 10
[D] Front	min. 10
[E] Rear	min. 10
Distance in cm (for operation in a tub)	
<b>Side</b>	
[A2] Top	can be located under a bench
[B] Left	min. 20
[C] Right	min. 20
[D] Front	min. 20
[E] Rear	min. 20

### 2.3.1 EMC-specific notes

**INFORMATION**

**Connecting cables in general**

Prerequisites for a failure-free operation of the temperature control units incl. their connections with external applications: Installation and wiring must be carried out professionally. Topics affected: "Electrical safety" and "EMC-compliant wiring".

**Cable lengths**

For flexible/fixe cable routing longer than 3 meters, the following must amongst other things be observed:

- Equipotential bonding, grounding (see also technical data sheet "Electromagnetic compatibility EMC")
- Compliance with "external" and/or "internal" lightning/overvoltage protection.
- Constructional protective measures, professional cable selection (UV resistance, steel pipe protection, etc.)

**Attention:**

The responsible body is responsible for compliance with national/international directives and laws. This also includes the testing of the installation/wiring required by law or standards.

These devices are suitable for the operation in "industrial electromagnetic environments". It meets the "immunity requirements" of the currently applicable EN61326-1, which are required for this environment.

It also meets the "interference emission requirements" for this environment. It is a **Group 1** and **Class A** unit according to the currently applicable EN55011.

**Group 1** specifies that high frequency (HF) is only used for the function of a device. **Class A** specifies the interference emission limits to be observed.

## 2.4 Installation conditions

 **WARNING**

Temperature control unit is connected to the power supply line

**DEATH FROM ELECTRICAL SHOCK BY DAMAGE TO THE POWER CABLE.**

- Do not put temperature control unit on power cable.

 **CAUTION**

Operating the temperature control unit fitted with castors without brakes activated

**CRUSHING OF LIMBS**

- Activate brakes on the wheels.

- Allow the temperature control unit to acclimate for about 2 hours when changing from a cold to a warm environment (or vice versa). Do not turn on the temperature control unit before!
- Install upright, stable and without tilt.
- Use a non-combustible, sealed subsurface.
- Keep environment clean: Prevent slip and trip hazards.
- Wheels must be locked after the installation, if installed!
- Spilled/leaked thermofluid must be disposed of immediately and properly. Follow the instructions on page 15, section »Proper disposal of resources and consumables« for the disposal of thermofluid and material.
- Observe the floor load bearing capacity for large units.
- Observe the ambient conditions.

## 2.5 Recommended temperature control and cooling water hoses



### Use of unsuitable/defective hoses and/or hose connections

#### INJURIES

- **Thermal fluid**
- Use appropriate hoses and/or hose connections.
- Check periodically for leaks and the quality of the hose and hose connections and take suitable measures (replace) as required.
- Isolate and protect temperature control hoses against contact/mechanical load.
- **Cooling water**
- Reinforced hoses must be used to satisfy tougher safety requirements.
- Shut off the cooling water supply to the temperature control unit even for shorter downtimes (e.g. overnight).



### Hot or cold thermal fluid and surfaces

#### BURNS TO LIMBS

- Avoid direct contact with the thermal fluids or the surfaces.
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).



### Uncontrolled formation of ice at the connections and hoses of the thermal fluid circuit

#### SLIP AND TRIP HAZARD

- If the temperature is controlled in the minus range, ice forms at the hoses and connections of the thermal fluid circuit. This occurs by condensing and freezing of atmospheric humidity.
- Check the strength of the ice formation. If too much ice is formed, this increases the risk of the temperature control unit tipping over. Secure the temperature control unit against tipping if this is the case.
- Check the ground below the ice formation for condensation water. Collect the condensation water with a suitable container or thoroughly remove it at regular intervals. You thus prevent the danger of slipping caused by condensation.

To connect applications, use only temperature control hoses that are compatible with the thermal fluid used. When selecting temperature control hoses, also pay attention to the temperature range in which the hoses are to be used.

- We recommend you use only temperature-insulated temperature control hoses with your temperature control unit. The responsible body is responsible for the insulation of connection valves.
- We **exclusively** recommend reinforced hoses for connecting to the cooling water supply. Cooling water and insulated temperature control hoses can be found in the Huber catalogue under Accessories.

## 2.6 Wrench sizes and torques

Note the wrench sizes that result for the pump connection on the temperature control unit. The following table lists the pump connections and the resulting wrench sizes, and torque values. A leak test must always be performed, and the connections tightened if necessary. The values of the maximum torque (see table) must **not** be exceeded.

Overview  
wrench sizes and  
torques

Connection	Sleeve nut wrench size	Connector wrench size	Recommended torques in Nm	Maximum torques in Nm
M16x1	19	17	20	24
M24x1.5	27	27	47	56
M30x1.5	36	32	79	93
	36	36	79	93
M38x1.5	46	46	130	153
G-thread (flat-sealing)	Adapt the torque to the material of the flat seal used. First hand-tighten the temperature control hose. When using adapters, do not overtighten the G-thread on the pump connection when connecting a temperature control hose. When connecting a temperature control hose to the adapter piece, secure the G thread against overwinding.			

## 2.7 Temperature control units with cooler

Valid heating circulators with >Cooling coil< [29]

### WARNING

Open electrical wires below the temperature control unit if the temperature falls below the dew point.

#### DEATH FROM ELECTRICAL SHOCK BY WATER ENTRY INTO THE ELECTRIC LINES.

- A temperature below the dew point may result in condensation in the temperature control unit and at the cooling water connections. The condensation is caused by high humidity at the cooling water-bearing components. The condensation exists the temperature control unit at the bottom.
- Electrical lines directly below the temperature control unit must be protected against liquid ingress.

### CAUTION

Use of unsuitable/defective hoses and/or hose connections

#### INJURIES

- **Thermal fluid**
- Use appropriate hoses and/or hose connections.
- Check periodically for leaks and the quality of the hose and hose connections and take suitable measures (replace) as required.
- Isolate and protect temperature control hoses against contact/mechanical load.
- **Cooling water**
- Reinforced hoses must be used to satisfy tougher safety requirements.
- Shut off the cooling water supply to the temperature control unit even for shorter downtimes (e.g. overnight).

### CAUTION

Steam or hot water escapes at the >cooling coil< [29]

#### SEVERE BURNS

- Do not use the >cooling coil< [29] at bath temperatures > 95 °C.
- Provide a greater flow rate through the >cooling coil< [29] at bath temperatures > 60 °C.

### NOTE

No protection against corrosion

#### DAMAGE TO THE TEMPERATURE CONTROL UNIT

- The addition of anti-corrosion agents is mandatory if salts (chlorides, bromide) have been added to the water circuit.
- Ensure that the materials used in the cooling water circuit are resistant with respect to the cooling water. See the data sheet from page 73 in section »Annex« for information on the materials used.
- Take suitable measures to maintain the warranty conditions.
- For information about water quality, see [www.huber-online.com](http://www.huber-online.com).

### NOTE

Usage of un-filtered river/sea or ocean water as cooling water

#### DAMAGE TO THE TEMPERATURE CONTROL UNIT

- Un-filtered river or sea water is not suitable for use as cooling water due to its contaminants.
- Use drinking water or filtered river or sea water for cooling.
- Sea water must not be used for water cooling.
- For information about water quality, see [www.huber-online.com](http://www.huber-online.com).

### INFORMATION

The >Cooling coil< [29] scales at cooling water temperatures > 60 °C according to the water quality. Descale regularly and on time. On page 15 observe Section »Proper disposal of resources and consumables«.

Observe all legal and other requirements specified by your water supplier when connecting to the public water system.

To reach temperatures close to the ambient temperature a >cooling coil< [29] is required. In some types of temperature control units the >cooling coil< [29] is pre-assembled and for other types it is available as an accessory (see the connection diagram on page 73 in section »Annex«).

## PROCEDURE

- Connect one of the connectors to your cooling water inlet and the other to your cooling water return flow (no preferred direction).
- Only use the >cooling coil< [29] at bath temperatures at less than 95 °C.
- Provide sufficient flow at high bath temperatures (between 60 °C and 95 °C) to prevent excessive




- scaling of the >cooling coil< [29] and the return line.
- Shut off the cooling water supply during heating. You will save energy and achieve shorter heating times.

## 2.8 Temperature control unit with batch

### 2.8.1 Operation as bath thermostat

Figure bath thermostat

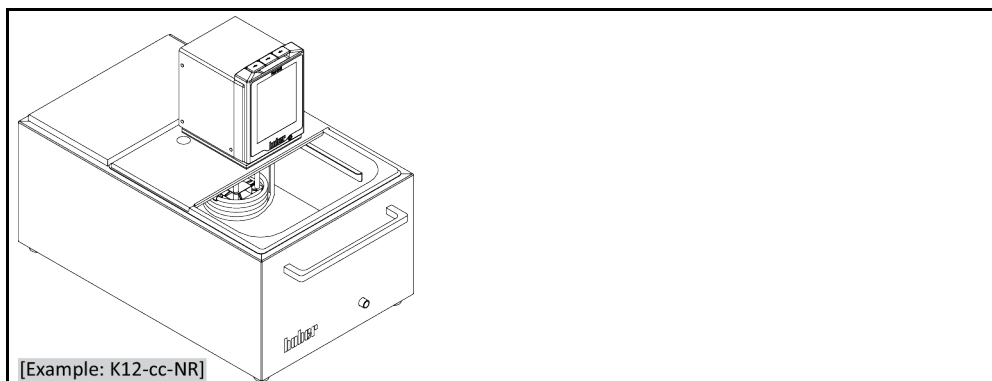

Note the volume displacement caused by a sample (e.g. Erlenmeyer flask). Place your sample into the empty bath. Only then fill in a sufficient amount of temperature control medium. Also note that the level of the temperature control medium drops when you remove the sample. This may cause a safety shutdown (low level protection) during an enabled thermoregulation. Therefore, switch off the temperature control unit beforehand.

## 2.9 Preparations for operation

### 2.9.1 Unscrewing/activating the leveling feet (if any)



**The leveling feet are not unscrewed/activated before switching on the temperature control unit  
DEATH OR SERIOUS INJURY DUE TO CRUSHING**

- The parking brakes must be activated at the casters (if any) and/or the leveling feet must be unscrewed/activated before the temperature control unit is put into operation.
- The temperature control unit may move if the parking brakes of the casters (if any) are not activated and/or the leveling feet are not unscrewed/activated.

Always unscrew/activate the leveling feet before switching on the temperature control unit. Uneven floors can be compensated by adjusting these leveling feet.

## PROCEDURE

- Verify that the parking brakes of the casters (if any) have been activated.
- Unscrew the leveling feet.
- Compensate uneven floors by adjusting these leveling feet, if necessary. Use a spirit level to horizontally align the temperature control unit.
- Tighten the lock screws on the leveling feet after aligning the temperature control unit. This prevents the leveling feet from changing their height during operation.

### 2.9.2 Externally closed and externally open applications

Using a pre-assembled pump adapter that is also available as an accessory, you can also control the temperature of an external application (e.g. reactor or open bath vessel). Externally open applications can run without interference only in conjunction with a DS level stabilizer (accessory). The DS level stabilizer compensates the differences in the pump (pressure capacity and throughput). If not already attached, please install the pump adapter. In an externally open application, please also install the DS level stabilizer on the externally open bath and observe on page 55 the section »Filling and venting the bath thermostat and the externally closed application« as well as the operation manual of the DS level stabilizer.

### 2.9.3 Opening/closing valves

Opening and closing valves



**INFORMATION**

**Opening valves:**  
Open valves by turning them counterclockwise (turn 90° left as far as it will go).  
**Close valves:**  
close valves by turning them clockwise (turn 90° right as far as it will go).

### PROCEDURE

- Check that all valves are closed.
- Close all valves by turning them clockwise (turn 90° right as far as it will go).

### 2.9.4 Installing collecting container

### PROCEDURE

- Install a suitable hose at the >overflow< [12] on the temperature control unit (if present). The hose must be compatible with the thermal fluid and the temperature.
- Put the other end of the hose in a suitable collecting container.

### 2.9.5 Connecting the functional earth

### PROCEDURE

- If required, connect the temperature control unit's >Functional earth terminal< [87] to the building's grounding point. Use a ground strap for this purpose. For the exact position and thread size please refer to the wiring diagram from page 73 in Section »Annex«.

## 2.10 Connecting externally closed application

The illustration "connection diagram" can be found on page 73 in section »Annex«.

### 2.10.1 Connecting an externally closed application


**NOTE**

**Overpressure in the application (e.g. > 0.5 bar (g) with glass apparatus)**

**MATERIAL DAMAGE TO THE APPLICATION**

- Provide an overpressure protective device to prevent damage to the application.
- Do not install valves/quick-release couplings in the feed/discharge lines from the temperature control unit to the application and from the application to the temperature control unit.
- **If valves/quick-release couplings are required:**
- Install burst disks on the application itself (at the feed and discharge lines).
- Install a bypass upstream of the valves/quick-release couplings for the application.
- Matching accessories (e.g. bypasses to reduce pressure) can be found in the Huber catalog.

Example: Connecting an externally closed application

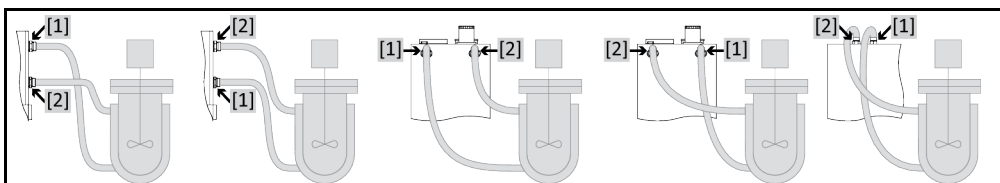

To enable your application to be operated correctly and eliminate air bubbles from the system, you must ensure that the >Circulation flow< [1] connection from the temperature control unit is at-

tached to the lower connection point of the application and the >Circulation return< [2] into the temperature control unit is attached to the higher connection point of the application.

## PROCEDURE

- Remove the screw plugs from the >Circulation flow< [1] and >Circulation return< [2] connections.
- Then connect your application to the temperature control unit using suitable thermal fluid hoses. The corresponding wrench sizes can be found in the table on page 23 in section »Wrench sizes and torques«.
- Check the connections for leaks.

### 2.11 Connecting to the power supply

#### INFORMATION

Based on local circumstances, it may be that you need to use an alternative power cable instead of the supplied original power cable. Do not use a power cable that is longer than **3 m** to be able to disconnect the temperature control unit at any time from the mains. Have the mains cable only replaced by a qualified electrician.

#### 2.11.1 Connection using socket with protective earth (PE)

##### DANGER

**Connecting to a power socket without protective earth (PE)**  
**MORTAL DANGER FROM ELECTRIC SHOCK**

- Always connect the temperature control unit to safety sockets (PE).

##### DANGER

**Damaged power cable/power cable connection**  
**MORTAL DANGER FROM ELECTRIC SHOCK**

- Do not start up the temperature control unit.
- Isolate the temperature control unit from the power supply.
- Have the power supply cable/power supply connection replaced and inspected by an electrician.
- Do not use a power cable that is longer than **3 m**.

#### NOTE

**Incorrect power supply connection**  
**DAMAGE TO THE TEMPERATURE CONTROL UNIT**

- Your building's existing power supply voltage and frequency must match the data provided on the rating plate of the temperature control unit.

#### INFORMATION

In case of uncertainties about an existing protective earth (PE), have the connection inspected by an electrician.

#### 2.11.2 Connection via hard wiring

##### DANGER

**Connection/adjustment to the power supply not carried out by an electrician**  
**MORTAL DANGER FROM ELECTRIC SHOCK**

- Have the connection/adjustment to the power supply carried out by an electrician.

##### DANGER

**Damaged power cable/power cable connection**  
**MORTAL DANGER FROM ELECTRIC SHOCK**

- Do not start up the temperature control unit.
- Isolate the temperature control unit from the power supply.
- Have the power supply cable/power supply connection replaced and inspected by an electrician.
- Do not use a power cable that is longer than **3 m**.

#### NOTE

**Incorrect power supply connection**  
**DAMAGE TO THE TEMPERATURE CONTROL UNIT**

- Your building's existing power supply voltage and frequency must match the data provided on the rating plate of the temperature control unit.

### 2.11.3 Converting the power supply connection

**INFORMATION**

To restore the factory settings, the connection values must be **re-entered**.

The following temperature control units can be retrofitted:  
CC-300BX, CC-315B(X)

- The temperature control units can be set to either 230 V/50 Hz or to 400 V/50 Hz.
- Install the temperature control unit at the location provided.
- The supply voltage must be selected during start-up for some temperature control units. This selection can be changed later in the Pilot ONE menu ([Category menu] > [Installation] > [Power supply]). The menu option "Voltage Supply" is visible only if the temperature control unit supports manually resetting the voltage supply. The supply voltage must be reselected if the temperature control unit has been reset to the factory settings.

## 3 Function description

### 3.1 Function description of the temperature control unit

#### 3.1.1 General functions

The **CC heating circulators** are **temperature control units** ideally suited for the thermoregulation of baths. This **controller generation** is characterized by the fact that there is only a **single hardware**. For a license fee you will receive an **activation key**, which allows you to quickly expand the functionality of e.g. a simple temperature control unit (basic version) to its **maximum functionality** (professional version) at any time.

Due to the powerful **heat technology**, **short heating times** can be achieved.

#### 3.1.2 Other functions

Due to the **speed-controlled pressure and suction pump**, the mixing of the thermal fluid can be optimally adapted to a given task.

The **self-optimizing cascade control** delivers **optimal controller results** for regulating levels up and down. Temperature control can be either a-periodic or with a slight overshoot (faster).

**Depending on the configuration level** both analog and digital interfaces, programmer and a Pt100 connection are available.

Your temperature control unit can be easily integrated into many laboratory automation systems. This is provided for by the existing **standard interfaces Ethernet, USB device and USB host at the "Pilot ONE"**. The optional Com.G@te allows you to expand your temperature control unit by **digital interfaces (RS232 and RS485), an analogue current loop interface (0/4 - 20 mA or 0 - 10 V) and various digital control options (In/Out)**.

The **removable control panel ("Pilot ONE")** can also be used as a **remote control**. Please contact your dealer or Huber Sales Department if you need an extension cable. The telephone number of the Huber Sales Department can be found on page 72 in Section »**Contact data**«.

A **connection jack for Pt100 process controller sensor** enables you to accomplish **external temperature control tasks** with ease (requires E-grade Exclusiv or Professional).

The temperature control unit can be retrofitted with an **integrated temperature ramp function** and an **internal temperature programmer** via an "E-grade" upgrade. The "E-grade" upgrade "**Exclusiv**" enables you to set and call 3 temperature control programs with 5 program steps each by using the integrated programmer. The "E-grade" upgrade "**Professional**" enables you to set and call 10 temperature control programs with a total of 100 program steps by using the integrated programmer.

The temperature control unit has **overtemperature protection to DIN EN 61010-2-010**, independent of the control circuit itself.

### 3.2 Information on the thermal fluids



**CAUTION**

**Non-compliance with the safety data sheet for the thermal fluid to be used**

#### **INJURIES**

- Risk of injury to the eyes, skin, respiratory tract.
- The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected.
- Observe the local regulations/work instructions.
- Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
- Danger of slipping because floor and work area are contaminated. Clean the work station and follow the instructions for the disposal of thermal fluid and material on page 15 in section »**Proper disposal of resources and consumables**«.

**NOTE**

**Non-compliance with the compatibility between the thermal fluid and your temperature control unit**

**MATERIAL DAMAGE**

- Observe the classification of your temperature control unit according to DIN 12876.
- Ensure the following materials are resistant with respect to the thermal fluid: Stainless steel 1.4301/ 1.4401 (V2A), copper, nickel, FKM, red bronze/brass, silver solder and plastic.
- The maximum viscosity of the thermal fluid must not exceed 50 mm<sup>2</sup>/s at the lowest working temperature!
- The maximum density of the thermal fluid may not exceed 1 kg/dm<sup>3</sup>!

**NOTE**

**Mixing different thermofluids in a thermal fluid circuit**

**PROPERTY DAMAGE**

- Do **not** mix different types of thermofluid (such as mineral oil, silicone oil, synthetic oil, water, etc.) in a thermofluid circuit.
- The thermal fluid circuit **must** be rinsed when changing from one type of thermal fluid to another. No residues of the previous type of thermal fluid may remain in the thermal fluid circuit.

Thermal fluid: Water

Designation	Specification
Calcium carbonate per liter	≤ 1.5 mmol/l; corresponds to a water hardness of: ≤ 8.4 °dH (soft)
PH value	between 6.0 and 8.5
Ultrapure water, distillates	Add 0.1 g of sodium carbonate (Na <sub>2</sub> CO <sub>3</sub> ) per liter
Not approved water	Distilled, deionized, demineralized, chloric, ferruginous, ammoniacal, or contaminated river water or sea water
Volume circulated (at least)	3 l/min.
<b>Thermal fluid: Water without ethylene glycol</b>	
Use	≥ 0 °C
<b>Thermal fluid: Water-ethylene glycol mixture</b>	
Use	< 0 °C
Thermal fluid composition	The mixture's temperature must be 10 K below the permissible min. temperature. For the permissible temperature range, refer to the datasheet from page 73 in Section »Annex«.

**INFORMATION**

For thermal fluids we recommend the media listed in the Huber catalog. The name of a thermal fluid is derived from its working temperature range and its viscosity at 25 °C.

### 3.3 To be noted when planning the test

**INFORMATION**

Also observe page 13 in section »Proper operation«.

The focus is on your application. Bear in mind that system performance is influenced by heat transfer, temperature, thermal fluid viscosity, volume flow, and flow speed.

- Make sure that the electrical connection is adequately dimensioned.
- The installation location of the temperature control unit should be selected so as to ensure adequate fresh air, even with water-cooled chillers.
- Select the thermal fluid to be used in such a way that it not only permits the minimum and maximum working temperature but is also suitable with regard to fire point, boiling point, and viscosity. In addition, the thermal fluid must be compatible with all the materials in your system.
- Basically, you should only use the thermal fluid recommended by the manufacturer and only within the usable temperature and pressure range.
- The application should be roughly at the same height of or below the temperature control unit if the thermoregulation is close to the boiling temperature of the thermal fluid.

- Fill the temperature control unit slowly, carefully and evenly. Wear the necessary personal protective equipment, such as goggles, heat-proof and chemical-resistant gloves, etc.

**INFORMATION**

For water-cooled temperature control units, please take the cooling water temperature necessary for perfect operation and the required differential pressure from the data sheet from page 73 onward in the Section »Annex«.

### 3.4 “Pilot ONE®” controller

The figure showing »“Pilot ONE”« can be found on page 6.  
The basic version of “Pilot ONE” (Basic) can be upgraded in three stages (from “Basic” to “Exclusive”, from “Exclusive” to “Professional”, and from “Professional” to “Explore”).

#### 3.4.1 Functional overview of “Pilot ONE®”

You can verify, and optionally upgrade, the delivery version of your temperature control unit using Pilot ONE, category “E-grade”.

Overview of the E-grade variants

Temperature control units / E-grade	E-grade Basic	E-grade Exclusive	E-grade Professional
Brewing thermostats	–	X	O
Unistat temperature control units	–	–	X
UniCAL	–	–	X
Other temperature control units	X	O	O
<b>E-grade “Explore”</b> (only for “Unistat”-series temperature control units) The E-grade includes E-grade “Professional” functionalities. In addition it includes: – Sensor calibration for return sensor – Return temperature display – Display of performance (estimated or calculated from volume flow) – Display mode “Explore” instead of “numerically large” – PB commands correspond to DV-E-grade, except: Temperatures only in 0.01 °C and volume flows only in 0.1 l/min.	–	–	O
<b>E-grade “DV-E-grade”</b> – All interface commands are enabled. – Temperatures in 0.001 C, volume flows in 0.001 l/min possible (cf. E-grade Explore) <b>Note:</b> Only the interface commands are enabled, not the corresponding menu items in Pilot ONE!	O	O	O
<b>E-grade “OPC-UA”</b> – OPC-UA interface via Ethernet. This E-grade additionally includes DV-E-grade functionalities.	O	O	O
X = Standard equipment, O = optional, – = not possible			

Overview of E-grade functions

Function	E-grade Basic	E-grade Exclusive	E-grade Professional
<b>Temperature control</b>			
Controller parameterization: predefined <sup>1</sup> / TAC <sup>2</sup>	X/–	–/X	–/X
Find parameters: Fast Identification / With Preliminary Test / Estimate Control Parameters	–/–/X	X/X/–	X/X/–

<sup>1</sup> TAC function available as a 30 day evaluation version.

<sup>2</sup> True Adaptive Control.

Function	E-grade Basic	E-grade Exclusive	E-grade Professional
Sensor calibration for control sensor <sup>1</sup> : x -point	2	5	5
Monitoring: Low level and overtemperature <sup>2</sup>	X	X	X
Adjustable alarm thresholds	X	X	X
VPC <sup>3</sup> (Variable Pressure Control)	X	X	X
Ventilation program	X	X	X
Automated compressor	X	X	X
Setpoint limitation	X	X	X
Program controller: x programs / max x steps	-/-	3/15	10/100
Ramp function: linear / non-linear	-/-	X/-	X/X
Temperature control mode: Internal and process	-	X	X
Adjustable max. heating and cooling capacity	-	X	X
<b>Display &amp; operation</b>			
Temperature display: 5.7" Touchscreen	X	X	X
Display mode: graphically / numerically large / Explore	-/X/-	-/X/-	X/X/-
Display resolution: 0.1 °C / 0.01 °C	X/-	X/X	X/X
Graphic display for temperature curves: Window, full screen and scalable	X	X	X
Calendar, date and time	X	X	X
Language: CZ, DE, EN, ES, FR, IT, JP, KO, PL, PT, RU, TR, ZH	X	X	X
Temperature format switchable: °C, °F and K	X	X	X
Display mode (screen) can be switched by swiping	X	X	X
Favorites menu	X	X	X
User menu (Administrator Level)	-	-	X
2. Setpoint	-	-	X
<b>Connections</b>			
Digital interface RS232	X	X	X
USB interfaces: Host und Device	X	X	X
Ethernet RJ45 interface	X	X	X
Pt100 external sensor connection	-	X	X
External control signal (ECS STANDBY <sup>4</sup> )	X	X	X
Programmable potential-free contact (ALARM <sup>5</sup> )	X	X	X
AIF (Analog Interface) 0/4–20 mA or 0–10 V <sup>6</sup>	X <sup>7</sup>	X	X
Digital interface RS485 <sup>1</sup>	X	X	X

<sup>1</sup> Internal Pt100 and External Pt100.

<sup>2</sup> For temperature control units with integrated overtemperature protection.

<sup>3</sup> For units with variable-speed pump or external bypass.

<sup>4</sup> Standard on Unistats, otherwise via optional Com.G@te or POKO/ECS interface.

<sup>5</sup> Standard on Unistats, otherwise via optional Com.G@te or POKO/ECS interface.

<sup>6</sup> Via optional Com.G@te.

<sup>7</sup> Restricted, see entries under "Temperature Control".



Function	E-grade Basic	E-grade Exclusive	E-grade Professional
<b>Comfort &amp; other</b>			
Visual / audible alarm signal	X	X	X
Autostart (Power Failure Auto)	X	X	X
Plug & Play technology	X	X	X
Technology glossary	X	X	X
Remote control / data visualization via Spy software	X	X	X
E-grade evaluation versions available (valid for 30 days)	X	X	X
Storage/loading temperature control programs onto USB stick	–	X	X
Copying settings	–	–	X
Service data recorders (flight recorder)	X	X	X
PB commands <sup>2</sup>	X	X	X
Communication Watchdog	–	–	X
Process data recording directly to a USB stick: Setpoint, Actual Value Internal and Actual Value Process / Heating Capacity %, Cooling Capacity %, and Pump Pressure / Pump Speed, and VPC Pressure	–/–/–	X/X/–	X/X/X
Calendar start	–	–	X

## 3.5 Clock/event function

### 3.5.1 Rechargeable accumulator

“Pilot ONE” is fitted with a clock that continues to run even when the temperature control unit is switched off. The energy required for this purpose is provided by a rechargeable accumulator, which is automatically charged when the temperature control unit is switched on. The accumulator is dimensioned so that the clock can also continue to run for prolonged switch-off intervals (up to several months). If, after extremely prolonged switch-off time, time and date have been deleted, leaving the temperature control unit switched on for a few hours will usually suffice (no temperature control required). During this time, you can reset the time and date.

If after switching the unit off and back on again, the previously set time and date re-appear, it can be safely assumed that the rechargeable accumulator is defective. In this case, contact Customer Support. The telephone number can be found on page 72 in section »Contact data«.

### 3.5.2 Programmable event function

The “Calendar Start” offers a programmable event function. This enables you to enter a time at which the event is repeatedly triggered on a daily basis (until the activity in the menu is reset). 2 event types are currently selectable:

#### 3.5.2.1 Event function “Alarm clock event”

Several acoustic signals are used.

#### 3.5.2.2 Event function “Program event”

After selecting “Program event” when configuring the event function, you will be prompted for the number of the program to be started. The program will be started automatically when the programmed event time is reached. If the temperature control unit is not active, this will also be started.

<sup>1</sup> Via optional Com.G@te.

<sup>2</sup> See Data Communication Manual. Anything controllable via the graphical user interface of the Pilot ONE can be controlled by PB commands.

### 3.6 Operation via the touch screen


The entire operation is via the >Touchscreen< [88]. These functions can be activated by tapping the displayed text boxes/icons once. This also changes the display.

**INFORMATION**

You can cancel the current dialog or dialog sequence at any time by pressing the "ESC" touch button. When canceling a dialog or dialog sequence, it may be necessary to confirm the cancellation again. When canceling a dialog sequence, settings made earlier in the dialog sequence are discarded. Check your already carried out settings and re-enter as needed.

### 3.7 Display instruments

Display instruments

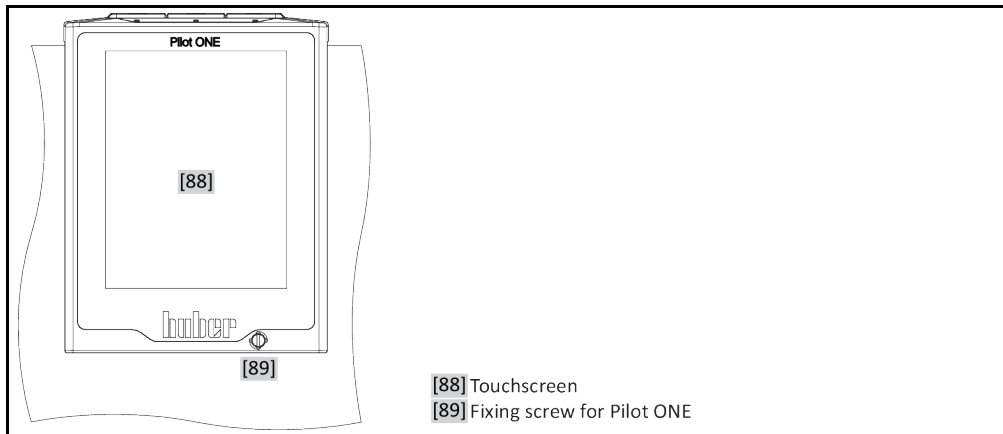

The following meters are available:


- >Touchscreen< [88]

#### 3.7.1 The touchscreen [88]

The most important display and operating instrument. Shows both standard variables (setpoint, actual value, setpoint thresholds...), and also menu guidance, error information output and operation.

### 3.8 Control instruments

Example "Control instruments"

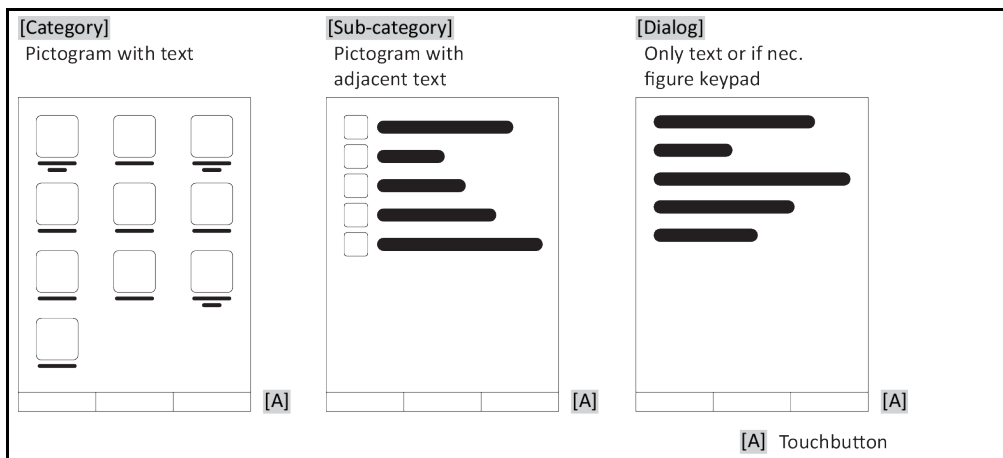

**INFORMATION**

To exit the “Categories Menu”, sub-categories, menu items, press the “Home” touch button (house) or the arrow. After 2 minutes of inactivity, the category/ sub-category or the Favourites menu is automatically closed and you return to the “Home” screen. Dialogs are **not** canceled/closed after 2 minutes of inactivity.

### 3.8.1 The touchbuttons

Depending on the situation, the touch buttons can be assigned different functions. For example:

- Select the “Home” screen (house)
  - Back (arrow to left)
  - Favourites (star)
  - Add to favorites (star with a plus sign)
  - Select the “Categories menu” (menu)
  - Confirm entry
  - Start/stop
- etc.

### 3.8.2 The categories

For clarity we have grouped the Operation and Setting of Pilot ONE in various categories. A category is selected by tapping it.

### 3.8.3 The sub-categories

The sub-categories are parts of a category. This is where you will find the entries that we have grouped together for you in the selected category. Not all the categories also contain sub-categories. Tap on a sub-category to select it.

### 3.8.4 The dialogs

Tapping on a category or sub-category displays its dialogs. Dialogues may appear e.g. as text, a numeric or a alphanumeric keyboard. Dialogs allow you for example to enter settings or start created temperature control programs. Within a dialog, selection must always be confirmed with the “OK” touch button. When cancelling a dialog with the “ESC” touch button it may be necessary to confirm the cancellation again.

## 3.9 Function examples


### 3.9.1 Display of software version

#### PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “System Settings”.
- Tap on the category “System Info”.
- Tap on the sub-category “Software Version”.

The software versions of the electronics will be displayed:

Display of software versions

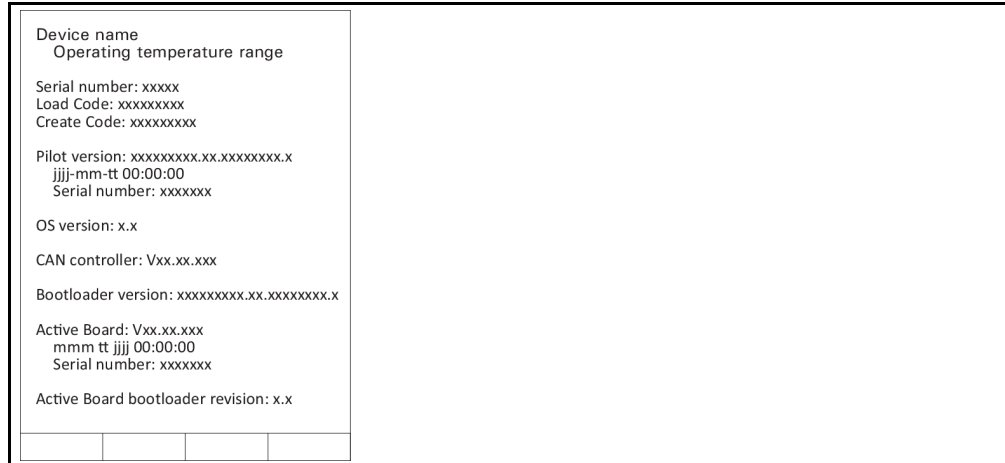

- Tap either the “ESC” or “OK” touchbutton. You will be returned to the previous screen.
- Tap on the “Home” touchbutton (house) to return to the “Home” screen.

### 3.9.2 Start & Stop

How to start and stop the temperature control process. Prerequisites: You have entered a setpoint.

## PROCEDURE

- Go to the “Home” screen.

#### Start

- Tap on the “Start” touchbutton.
- Confirm the start of temperature control by tapping on “OK”.  
The correct selection will be displayed graphically and temperature control will start immediately. If tapping on “OK” is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the “Home” screen again. Try to start the temperature control unit again.

#### Stop

- Tap on the “Stop” touchbutton.
- Confirm the stop of temperature control by tapping on “OK”.  
The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on “OK” is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the “Home” screen again. Try to stop the temperature control unit again.

#### INFORMATION

With the “Stop” touchbutton, you can also stop Degassing, Venting and Circulation. The requirement for this is that the corresponding task is active.

### 3.9.3 Copying the settings to a data carrier

#### Only valid for E-grade Professional

The current settings of the temperature control unit are saved as a file to a connected USB flash drive. Use this file to restore the settings of a temperature control unit or to copy them to another temperature control unit. The models of the temperature control units used must be identical to copy the settings between temperature control units. You cannot transfer data between different models.

You can also transfer the set value of the overtemperature protection to another temperature control unit. You must **always** verify this value on the respective temperature control unit and adjust it if necessary. Make sure that the value of the overtemperature protection has been adapted to the thermal fluid used.

**3.9.3.1 Saving to a USB flash drive**

**PROCEDURE**

- Insert a USB flash drive into the “USB 2.0 interface host”. The USB flash drive must have at least 1 MB of free space.
- Go to the “Categories Menu”.
- Tap on the category “System Settings”.
- Tap on the category “Copy Settings”.
- Tap on the dialog entry “Save to USB flash drive”.
- Confirm your choice by tapping on “OK”. The content of the USB flash drive is displayed. If necessary, select the location (folder) on the USB flash drive.
- Confirm the selection of the memory location by tapping on “OK”.
- Enter a name for the file to be saved. Confirm the entry by tapping on “OK”. Alternatively, you can accept the proposed name by tapping on “OK”.
- Read the message and confirm by tapping on “OK”. The file with the settings was created on the USB flash drive.
- Remove the USB flash drive from the temperature control unit.

**3.9.3.2 Loading from a USB flash drive**

**PROCEDURE**

- Insert a USB flash drive with the saved file into the “USB 2.0 interface host”.
- Go to the “Categories Menu”.
- Tap on the category “System Settings”.
- Tap on the category “Copy Settings”.
- Tap on the dialog entry “Load from USB flash drive”.
- Confirm your choice by tapping on “OK”. The content of the USB flash drive is displayed.
- Select the file you want to load.
- Confirm the selection of the file by tapping on “OK”.
- Select the setting group to be loaded from the list. A multiple selection is possible.
- Confirm your choice by tapping on “OK”.
- Read the message on the >Touchscreen< [88].
- Switch off the temperature control unit. The settings have been loaded to the temperature control unit.
- Remove the USB flash drive from the temperature control unit.

**3.9.4 Restore factory settings**

Use this function to reset the temperature control unit to various basic states. This is useful if you want to reverse various settings relatively quickly.

**INFORMATION**

Restoring the factory settings is **only** possible if the temperature control unit is **not** carrying out a task. If a task is active, do not switch off the temperature control unit until the application allows this. Restoring the factory settings cannot be revoked. Depending on the type of resetting to factory settings performed you may have to enter parameters (process safety, thermal fluid used, cut-out value etc.) again.

“X” = Value is reset, “-” = Value is not reset

(A) = All together; (B) = Device parameters without OT; (C) = Device Parameter; (D) Program Controller; (E) = Menu; (F) = Com.G@te

Description	(A)	(B)	(C)	(D)	(E)	(F)
<b>Acoustics</b>	X	X	X	-	-	-
<b>Display functions</b> Warnings; Messages; Inactive Menu Items; Resolution Temperature	X	X	X	-	-	-
<b>Display modes</b>	X	X	X	-	-	-
<b>AutoStart response</b>	X	X	X	-	-	-
<b>Limitations</b> DeltaT limit; max. heating variable; max. cooling variable; max. permissible current consumption (for 230 V temperature control units)	X	X	X	-	-	-

Description	(A)	(B)	(C)	(D)	(E)	(F)
<b>Com.G@te (only if Com.G@te is connected)</b> Analog Interface (configuration input)						
Meaning of AIF input; response on cable break; measurement range setting	X	X	X	-	-	X
Calibration	-	-	-	-	-	X
In case of analog error; current / voltage switching; filter constant	X	X	X	-	-	X
<b>Configuration Output</b> AIF output variable; Output range setting	X	X	X	-	-	X
Calibration	-	-	-	-	-	X
RS232/RS485 (Hardware RS, Baud rate, device address); ECS Standby; POKO Alarm	X	X	X	-	-	-
<b>Properties (other)</b> Change thermal fluid (thermal fluid, liter indicator, bypass usage); cooling bath / cooling (CC-E only)	X	X	X	-	-	-
<b>Sensor calibration</b>	-	-	-	-	-	-
<b>Automated compressor</b>	X	X	X	-	-	-
<b>Enter program</b>	X	-	-	X	-	-
<b>Pump settings</b> Pump speed setpoint; Pump pressure setpoint; Control Mode	X	X	X	-	-	-
<b>Fluid properties</b> Thermal fluid; liter indicator	X	X	X	-	-	-
Configuration manual parameters						
<b>Protective function</b> High Limit Alarm Internal; Low Limit Alarm Internal; High Limit Alarm Process; Low Limit Alarm Process; Hydrostatic Correction; Warning Time Limit (CC-E only); Min. Level (for temperature control units with analogue level sensor); Max. Level (for temperature controllers with analogue level sensor); Freeze Protection (optional)	X	X	X	-	-	-
<b>Setpoint</b>	X	X	X	-	-	-
<b>Setpoint limitation</b> Min. Setpoint; Max. Setpoint	X	X	X	-	-	-
<b>Language</b>	X	X	X	-	-	-
<b>Temperature format</b>	X	X	X	-	-	-
<b>Temperature control mode</b>	X	X	X	-	-	-
<b>Favourites menu</b>	X	-	-	-	X	-
<b>Pilot ONE Ethernet</b> IP address; Subnet Mask; Remote Access	X	X	X	-	-	-
<b>Time format</b>	X	X	X	-	-	-
<b>2. Setpoint</b>	X	X	X	-	-	-

### 3.9.4.1 Restore to factory settings without overtemperature protection

## PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "System Settings".
- Tap on the category "Factory Settings".
- Choose between the dialog entries "Unit Control Data without OT", "Menu", "Programmer" and "Com.G@te". These entries do **not** restore the overtemperature protection. Tap on the required dialog entry.

- Tap on the “OK” touchbutton to confirm your choice.
- Read the message displayed. Tapping on “Yes” restores the factory settings, tapping on “No” cancels the procedure. The message “Restart system!” appears on the >Touchscreen< [88].
- Switch the temperature control unit off. The selected control data have been reset.

### 3.9.4.2 Restore to factory settings including overtemperature protection

## PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “System Settings”.
- Tap on the category “Factory Settings”.
- Choose between the dialog entries “Unit Control Data” and “All together”. These entries **also** reset the overtemperature protection. Tap on the required dialog entry.
- Tap on the “OK” touchbutton to confirm your choice.
- Read the message displayed. Tapping on “Yes” restores the factory settings, tapping on “No” cancels the procedure.

## INFORMATION

In the following dialog, enter the overtemperature protection suitable to the thermal fluid used. To reset the cut-out values of the overtemperature protection to the factory setting from within the controller of Pilot ONE, enter for “Heater” 35 °C and for the expansion vessel 45 °C. “Process Safety” is factory set to “Stop” and is automatically reset to “Stop” when restoring the factory settings.

- Enter the temperature unit you want to use with Pilot ONE. The available choices are “Celsius (°C)”, “Kelvin (K)” and “Fahrenheit (°F)”.
- Confirm your choice by tapping on “OK”.
- Read the message and confirm by tapping on “OK”.
- Read the safety warning and confirm by tapping on “OK”.
- Read the Note and confirm by tapping on “OK”.
- Tap on the temperature format set by you in the controller (green text).
- Enter the displayed red safety code using the numeric keypad that appears.
- Enter the value 35 °C for the “Heating OT value” using the numeric keypad that appears.
- Confirm your entry by tapping on “OK”.
- Read the message and confirm by tapping on “OK”.
- Read the safety warning and confirm by tapping on “OK”.
- Tap on the temperature format set by you in the controller (green text).
- Enter the displayed red safety code using the numeric keypad that appears.
- Enter the value 45 °C for the “Exp. vessel OT value” using the numeric keypad that appears.
- Confirm your entry by tapping on “OK”. The message “Restart system!” appears on the >Touchscreen< [88].
- Switch off the temperature control unit. The temperature control unit is reset.

## 4 Setup mode

### 4.1 Setup mode

#### CAUTION

**Moving the temperature control unit during operation**

**SERIOUS BURNS/FREEZING OF THE HOUSING PARTS/ESCAPING THERMAL FLUID**

➤ Do not move temperature control units that are in operation.

#### 4.1.1 Turning on the temperature control unit

##### PROCEDURE

- Switch on the temperature control unit using the **>mains switch<** [37].  
First, a system test is conducted to test the full functionality of the temperature control unit. All sensors, the all pole breaking power disconnect relay for the power supply and the power electronics of the main heater and the main heater itself are checked. A message appears on the **>touchscreen<** [88] in case of an error or if a warning is issued. If in doubt contact Customer Support (the phone number can be found on page 72 in section **»Contact data«**).

#### INFORMATION

The following data entry is required only for:

- a) initial startup
- b) resetting the temperature control unit to the factory settings  
(see page 37, section **»Restore factory settings«**).

- After switching on the temperature control unit, tap on the required system language.
- Confirm your choice by tapping on "OK".
- Tap on the thermal fluid used.
- Confirm your choice by tapping on "OK".
- Read the message and confirm by tapping on "OK".
- Enter the volume using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Tap on the cooling bath used. This selection list is displayed only when the Pilot ONE was connected to a cooling bath.
- Confirm your choice by tapping on "OK".

#### INFORMATION

If your temperature control unit is not integrated into a network, confirm the preset IP address (0.0.0.0) by tapping on "OK". This skips the network settings.

- Enter the required IP address using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the required Subnet Mask using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Tap on the required remote control mode.
- Confirm your choice by tapping on "OK".

#### 4.1.2 Turning off the temperature control unit

##### PROCEDURE

- Warm the thermofluid to room temperature.
- Stop the thermoregulation.
- Switch off the temperature control unit using the **>Mains switch<** [37].



### 4.1.3 Setting the overtemperature protection



**The overtemperature protection is not correctly set to the thermal fluid used**

#### **MORTAL DANGER FROM FIRE**

- The overtemperature protection must be correctly set to the thermal fluid you are using.
- Always observe the safety data sheet of the thermal fluid.
- Different working temperature ranges may apply for open and closed systems.
- **For Huber thermal fluids:**
- For Huber thermal fluids, the maximum usable working temperature range is already indicated.
- Set the cut-out value of the overtemperature protection to the **upper temperature limit of the working range** of the thermal fluid.
- **For thermal fluids of other manufacturers:**
- Set the cut-out value of the overtemperature protection at least **25 K below the fire point** of the thermal fluid.

#### **INFORMATION**

The maximum settable cut-out value of the OT corresponds, for Huber thermal fluids, to the specified upper working temperature of the thermal fluid. The usable working temperature range may be less, if the overtemperature protection is correctly set. When controlling the temperature at the upper working temperature limit, the OT may be triggered due to tolerances.

#### 4.1.3.1 General information on the overtemperature protection

The overtemperature protection is a device in the temperature control unit that works independently of the controller. The software and hardware is designed to test key functions and operating states during a self-test after the power supply has been switched on. If faults are detected, the enabling of the electrical assemblies in the temperature control unit is blocked. During operation, the sensors are tested for short-circuit and interruption.

The bath or flow temperature is monitored to ensure the safety of your system. It is set immediately after you have filled the system with thermal fluid.

Our temperature control units not only offer the possibility to specify the cut-out value of the overtemperature protection but also provide the possibility to specify the shutdown mode of the temperature control unit. In a typical setting the temperature control unit switches off both the temperature control as well as the circulation (**stop according to DIN EN 61010**) after reaching the cut-out value. This monitors a possible defect in the control of the heater. A strong exothermic reaction near the cut-out value can also, under certain circumstances, result in a shutdown of the temperature control unit. In this case, however, a shutdown would be fatal.

Our temperature control units provide the possibility to work with the shutdown mode **Process Safety**. In this mode, temperature control (cooling) and circulation work as always. This is therefore a possibility to respond to exothermic reactions.

#### **INFORMATION**

The default shutdown mode of the overtemperature protection is **“Stop according to DIN EN 61010”**. When resetting to factory settings, the overtemperature protection is **reset** to the default shutdown mode **“Stop according to DIN EN 61010”**!

Upon delivery, the cut-out value of the overtemperature protection is set to 35 °C. An alarm is triggered by the temperature control unit shortly after turning on the power if the temperature of the thermal fluid just filled is higher than the cut-out value set for the overtemperature protection. Set the overtemperature protection to the thermal fluid you are using.

When setting a new cut-out value for the overtemperature protection, you will be prompted to enter a randomly generated and displayed code via the numeric keypad displayed. Only after successful entry will you be able to change the cut-out value.

#### 4.1.3.2 Setting “OT limit: heating”

### **PROCEDURE**

- Go to the “Categories Menu”.
- Tap on the category “Safety”.
- Tap on the category “Overtemperature”.

- Read the safety warning and confirm by tapping on “OK”.
- Read the Note and confirm by tapping on “OK”.
- Tap on the dialog entry “OT limit: heating”.
- Confirm your choice by tapping on “OK”.
- Tap on the temperature format set by you in the controller (green text).
- Enter the displayed red safety code using the number keypad that appears.
- Enter the value for the “Heating OT value” using the numeric keypad that appears. This value must be 25 K below the fire point of the thermal fluid you are using.
- Confirm your entry by tapping on “OK”.

#### 4.1.3.3 Setting “Process Safety”

You have two options.

##### “Stop” in accordance with DIN EN 61010

Once the overtemperature protection cut-out value is reached, the temperature control unit (heating, cooling circuit and circulation pump) switches off (default setting).

##### “Process Safety”

Once the overtemperature protection cut-out value is reached, the heater switches off and the cooling circuit and circulation pump continue to operate. In case of an emergency (possible exothermic reaction), the full cooling capacity is thus available. Please make sure that the automated compressor is set to **Permanently ON** ([System Settings] > [Power/ECO Settings] > [Compressor ON/OFF/AUTO] > [Permanently ON]).

## PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “Safety”.
- Tap on the category “Process Safety”.
- Read the safety warning and confirm by tapping on “OK”.
- Choose between the modes “Stop” and “process safety”.
- Confirm your choice by tapping on “OK”.

#### 4.1.3.4 Monitoring via “Display OT values”

## PROCEDURE

- Go to the “Categories Menu”.
  - Tap on the category “Safety”.
  - Tap on the category “Display OT values”.
- You will receive an overview of the current measured temperature value of the overtemperature protection sensor, the set cut-out values and the set cut-out mode (Process Safety). Some temperature control units have 2 overtemperature protection sensors, and so 2 values are shown for these temperature control units.
- Tap on the “OK” touchbutton after you have read/checked the information.

#### 4.1.4 Testing overtemperature protection for functionality



**Overtemperature protection (OT) does not trip**

##### **MORTAL DANGER FROM FIRE**

- Test the response of the device every month and after each change of the thermal fluid in order to assure proper functioning.

##### **NOTE**

**The steps below are carried out without permanent monitoring of the temperature control unit DAMAGE TO AND IN THE VICINITY OF THE TEMPERATURE CONTROL UNIT**

- The following actions may only be carried out while constantly monitoring the temperature control unit and the application!

**INFORMATION**

Run the test only if the temperature of the thermal fluid used is about 20 °C. You must **NOT leave the temperature control unit unattended** during the entire test of the overtemperature protection.

The Pilot ONE controller contains a description for carrying out the overtemperature protection test.

**PROCEDURE**

- Go to the “Categories Menu”.
- Tap on the category “Safety”.
- Tap on the category “OT Test”.
- The description for carrying out the overtemperature protection test is displayed.
- Tap on the “OK” touch button after you have read the information.

**4.1.5 Adjusting the Delta T limiter****NOTE**

**The Delta T limiter has not been adjusted to the used glass apparatus**

**DAMAGE CAUSED BY GLASS APPARATUS BURSTING**

- Adjust the Delta T value to your application.

**INFORMATION**

The temperature dynamics inside the reactor/process temperature are dictated by the flow temperature. A differential temperature (Delta T) occurs between the flow temperature and the temperature inside the reactor. The greater the permitted Delta T, the better the energy transmission and hence the faster the speed reaches the setpoint. However, damage could result if the temperature difference limits are exceeded (bursting of the application e.g. glass apparatus). This difference in temperature may have to be restricted depending on the application (glass apparatus).

**4.1.5.1 Changing the Delta T limiter****PROCEDURE**

- Go to the “Categories Menu”.
- Tap on the category “Protection Options”.
- Tap on the sub-category “Delta T limiter”.
- Set the value of Delta T in line with the glass apparatus.
- Confirm your entry by tapping on the “OK” touchbutton.

**4.2 The temperature control circuit**


Every temperature control unit with a Pilot ONE control panel has its own PID controller for internal and process temperature control. For many temperature control tasks it is sufficient to use the factory-set control parameter. Our years of experience and current developments in control technology are used in these control parameters.

If a process control system is used, ideally send the specification of the temperature setpoint digitally to the temperature control unit. For this purpose, an Ethernet and USB port is available at the Pilot ONE and a RS232 interface at the temperature control unit. With the optional Com.G@te, an additional RS485 interface is added to your temperature control unit. Optionally you can integrate the temperature control unit in a Profibus environment. For more information, see the pages from 60, section »**Interfaces and software update**«.

**INFORMATION**

The capacity adjustment of the temperature control unit is optimized so that the specified processes are run through in the fastest possible time. This increases the productivity of the system as a whole and thus saves energy in the long term.

Diagram of temperature control circuit

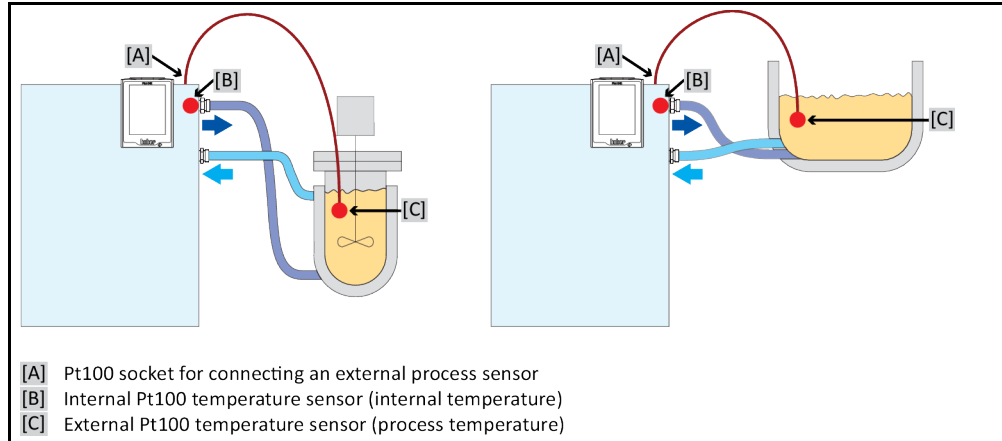

### 4.2.1 Select temperature control: Internal or process


#### PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “Temperature Control”.
- Tap on the category “Process/Internal”.
- Choose between the entries “Internal” and “Process (Cascade)”.
- Confirm your choice by tapping on “OK”.

### 4.2.2 Temperature control to internal temperature

With internal temperature control, a control circuit is used to control the temperature at the internal Pt100 temperature sensor. This Pt100 temperature sensor is built into the device and is located close to thermal fluid outlet (forward flow) or in the bath tank.

Representation of an optimum internal temperature control

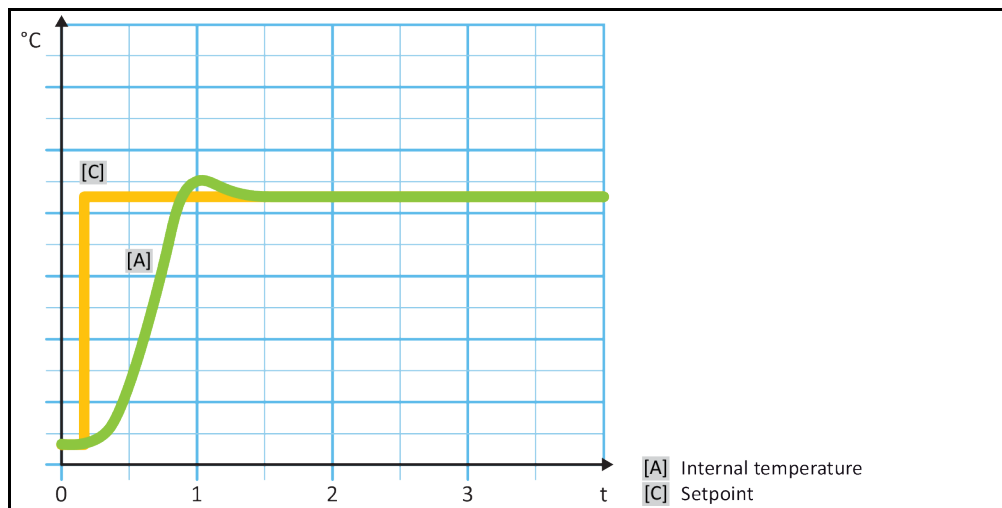

### 4.2.3 Temperature control to process temperature


Certain temperature control tasks require that the temperature is recorded elsewhere for best results than described. Setting the temperature to process temperature makes alternatives available. When setting the temperature to process temperature, an additionally connected external Pt100 temperature sensor is used in conjunction with a master controller (cascade controller). The internal sensor at the supply line is integrated with the slave controller. This temperature control method is used e. g. for the thermostatic control of jacket vessels. The setpoint setting is valid for the process controller. It in turn calculates a target value for the internal controller to optimally adjust the process setpoint.

**NOTE**

**Incorrect installation of the process sensor (Pt100)**  
**MATERIAL DAMAGE CAUSED BY FAULTY TEMPERATURE CONTROL**

- Disruption of measured value recording due to static build-up.
- The process sensor (Pt100) must have a screened supply cable.
- If the sensor tube is metallic, take care to avoid ground loops.
- The connection cable should not be unnecessarily long.
- Make sure that the process sensor is properly attached at the measurement point and that there is good thermal coupling
- The sensor itself must have good insulation from the screen or the protective grounding ( $R > 20 \text{ M}\Omega$ ).

Representation of an optimum process temperature control

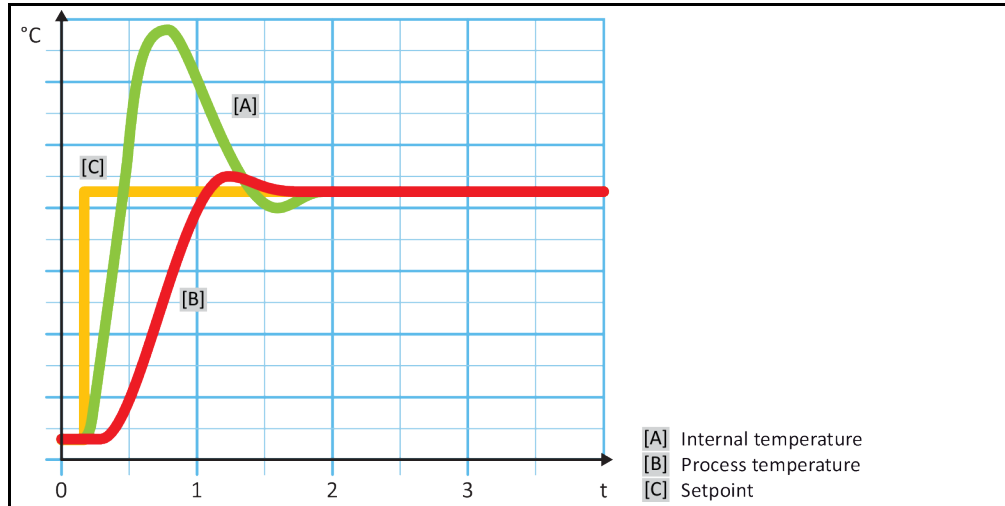

#### 4.2.4 Delta T limiter

The Delta T limiter is a part of the temperature control which protects the system or the process. The Delta T limiter is given a limit value. The Delta T limiter then reacts when the limit value is reached during heating or cooling.

The temperature control mode “Process (Cascade)” evaluates the temperature difference between the flow temperature and the process temperature. The default setting of the limit value is 100 K. If the limit value and the temperature sensor are set properly, the load limits, e.g. from a glass apparatus, are not exceeded. As the limit value is approached, the cooling or heating capacity is adjusted. The DeltaT limiter is **not** a safety device.

#### 4.2.5 Monitoring the Pt100 temperature sensors

The Pt100 temperature sensors are constantly monitored for their electrical status. If the status “Sensor faulty” occurs during temperature control, the temperature control process is stopped immediately and a device message is displayed. This applies for all temperature sensors connected within the temperature control unit.

#### 4.2.6 Optimum control parameters for optimum temperature control

If the adjustment of the temperature does not correspond to the quality of the illustrations shown above, you can adjust the control parameters. With Huber temperature control units, there are various ways of finding the optimum control parameters. Depending on the facilities of the temperature control unit, you can choose the following processes:

- Use ex-factory parameters (standard)
- Estimate Control Parameters (only useful for bath thermostats with E-grade Basic and internal control)
- Fast Identification (from E-grade Exclusive)
- With Preliminary Test (from E-grade Exclusive)

## 4.2.7 Sub-category: "Select auto/expert mode".

**NOTE**

Use of the "Expert mode" without a thorough knowledge of I&C technology.

**MATERIAL DAMAGE TO THE APPLICATION**

- Only use this mode if you have a thorough knowledge of I&C technology.

Here you can select whether the control parameters are set in the "Automatic mode" or in the "Expert mode". To enter settings in "Expert mode" requires in-depth knowledge of process measuring and control technology. Incorrect or insufficient settings can severely impair the function of the temperature control unit.

**INFORMATION**

In "Expert mode", "Configuration auto" is deactivated and only "Configuration manual" is possible.

This is how to change mode:

## PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Select auto/expert".
- Choose between the dialog entries "Automatic mode" and "Expert mode".
- Confirm your choice by tapping on "OK".

## 4.2.8 Sub-category: "Configuration auto"

**INFORMATION**

You can only select this menu entry if "Automatic mode" is set.

### 4.2.8.1 Sub-category: "Find parameters"

#### 4.2.8.1.1 Dialog entry: "Fast Identification"

Considering the small effort, the "Fast Identification" function of the controlled system provides you with very quick and reliable adapted control parameters. These control parameters achieve a fast and very accurate tuning performance. The more complex but also more precise identification "With Preliminary Test" is only very rarely required.

**INFORMATION**

Do not make any changes to the temperature control unit and to the application (temperature control unit / external application) after your system has been started. Changes include e.g. filling / draining process chamber, changing the mixer speed, change of position of Pt100 process control sensor etc.

## PROCEDURE

- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- Do not make any changes to the temperature control unit and the application during parameter determination.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Find parameters".
- Tap on the dialog entry "Fast Identification".
- Confirm your choice by tapping on "OK".
- Read the message displayed and confirm it by tapping on "OK".
- Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK".
- Choose between the dialog entries "Internal" and "Process (Cascade)".
- Confirm your choice by tapping on "OK".

- Enter a new setpoint using the number keypad that appears. This should be at least 10 K away from the current setpoint.
- Confirm your entry by tapping on "OK". Setting of the control parameters using "Fast Identification" starts and after a while a message appears on the display.
- Read the message displayed and confirm it by tapping on "OK".

#### 4.2.8.1.2 Dialog entry: "With Preliminary Test"

In some complex applications, the "Fast Identification" of the control circuit may not yet lead to optimum control. This can happen, in particular, if the hydraulic build-up cannot be adjusted to the necessary circulation volumes. Further information on this is provided on page 30 in the Section »**To be noted when planning the test**«.

A further optimization of control behavior can be achieved if you select the control parameterization "With Preliminary Test". With this, the control parameters within the set limits of the minimum and maximum setpoint are determined. Temperature control is then also carried out to the setpoint limits in some circumstances.

### INFORMATION

Before starting the automatic controller parameterization, verify the correct setting of the minimum and maximum setpoint. A limitation to the actually used operating temperature range is advantageous. Do not make any changes to the temperature control unit and to the application (temperature control unit / external application) after your system has been started. Changes include e.g. filling / draining process chamber, changing the mixerspeed, change of position of Pt100 process sensor etc.

Since the working temperature range can sometimes be very large, finding the parameters takes correspondingly longer in this mode. The control defines up to three temperature setpoints and automatically processes them one after another. One of these is below room temperature, one roughly at room temperature and one above room temperature, if the setpoint limits allow this.

## PROCEDURE

- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- Do not make any changes to the temperature control unit and the application during parameter determination.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Find parameters".
- Tap on the dialog entry "With Preliminary Test".
- Confirm your choice by tapping on "OK".
- Read the message displayed and confirm it by tapping on "OK".
- Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK".
- Choose between the dialog entries "Internal" and "Process (Cascade)".
- Confirm your choice by tapping on "OK". Setting the control parameters using "With Preliminary Test" starts and a message appears on the touchscreen after a while.
- Read the message displayed and confirm it by tapping on "OK".

#### 4.2.8.1.3 Dialog entry: "Estimate Control Parameters"

Even with the simple temperature control units, we offer a further advantage in comparison with comparable bath thermostats available on the market. You can modify an existing control parameter by entering the thermal fluid used and the quantity of thermal fluid. This version is available with bath thermostats without connected external application.

### INFORMATION

All the necessary technical data for the thermal fluids listed for Pilot ONE are stored in the controller. If the thermal fluid that you use is not shown in the list, choose the most similar thermal fluid with regard to temperature range and viscosity. Further information about this is provided on page 29 in the Section »**Information on the thermal fluids**«.

## PROCEDURE

- Before setting the control parameters, make sure that the temperature control unit has reached the set setpoint and has been controlling the temperature at this setpoint for a few minutes. Do not stop the temperature control.
- Do not make any changes to the temperature control unit and the application during parameter determination.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Find parameters".
- Tap on the dialog entry "Estimate Control Parameters".
- Confirm your choice by tapping on "OK".
- Read the message and confirm it by tapping on "Yes". The current control parameters will be overwritten.
- Select the thermal fluid used from the list.
- Confirm your choice by tapping on "OK".
- Read the message displayed and confirm it by tapping on "OK".
- Enter the filling volume using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Tap twice on the "Arrow" touchbutton to return to the category "Temperature Control".
- Tap on the category "Process/Internal".
- Choose between the dialog entries "Internal" and "Process (Cascade)".
- Tap on "OK" to confirm your choice.

### 4.2.8.2 Sub-category: "Control Dynamics"


You can choose between a faster tuning performance with a possible and accepted small temperature overshoot and a tuning performance without temperature overshoot. The default setting is "**Fast, small overshoot**".

The overshoot always refers to the leading temperature. Example: The process temperature, if activated, is the leading temperature. In contrast, the bath or flow temperature must always be leading the process temperature. In order to achieve the best possible energy transfer, the greatest possible temperature difference between the bath or flow temperature and the process temperature is necessary (see figure "»Representation of an optimum process temperature control«" from page 44 in section »Temperature control to process temperature«). This always requires the maximum possible thermal fluid flow rate. In the setting "**Faster, small overshoot**", the combination of high thermal fluid flow rate and the superbly designed control electronics hardly ever results in an overshoot of the process temperature, while reaching the setpoint as quickly as possible. The opposite to the mode "**Faster, small overshoot**" is available in the setting "**Without overshoot**". The target temperature is approached more cautiously and thus aperiodic. It takes longer to adjust to the selected setpoint. The statement "without overshooting" is only valid with little external disturbance influence. Please note the specifications on page 30 in section »To be noted when planning the test«.

#### INFORMATION


The adjustment behavior can be selected at any time without reactivating the "Find control parameters".

Internal, aperiodic temperature control

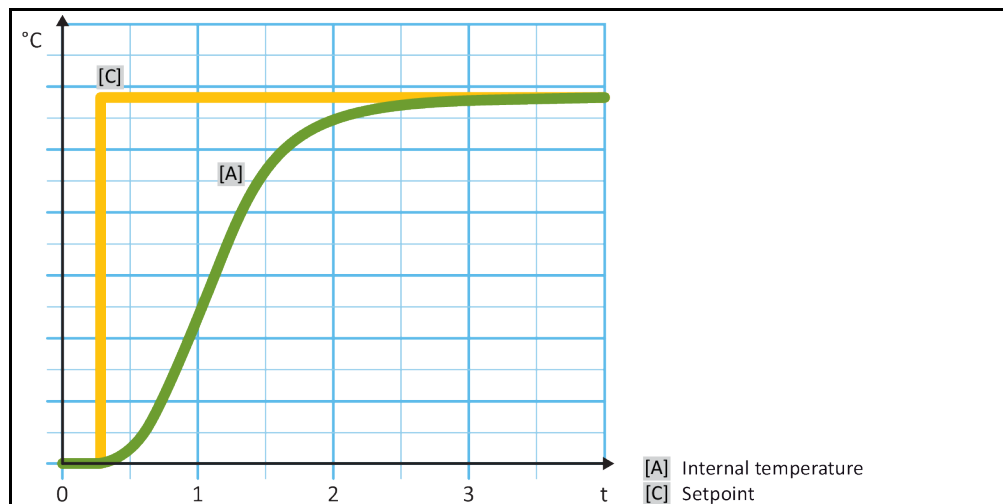



Internal, dynamic temperature control with possible overshooting of the temperature

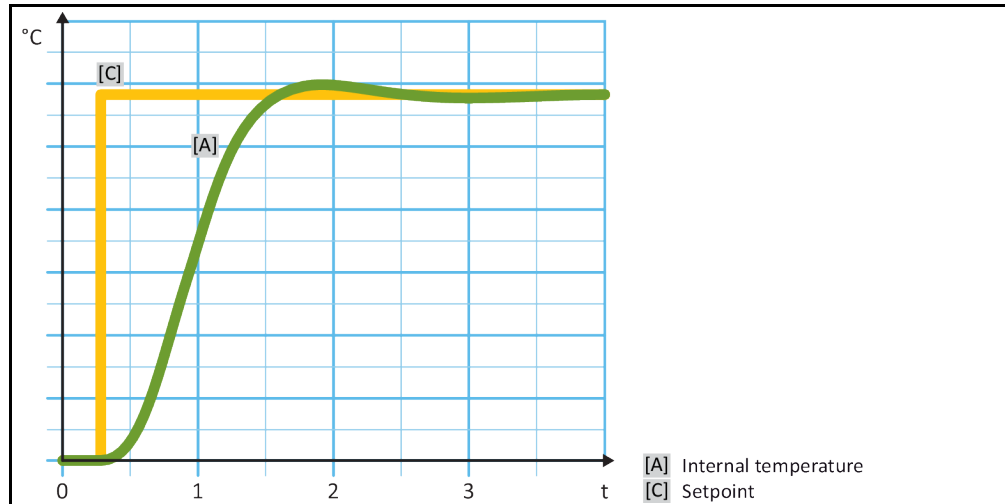

## PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “Temperature Control”.
- Tap on the category “TAC/Manual”.
- Tap on the sub-category “Configuration auto”.
- Tap on the sub-category “Control Dynamics”.
- Choose between the dialog entries “Fast, small overshoot” and “Without overshoot”.
- Confirm your choice by tapping on “OK”.

### 4.2.8.3 Sub-category: “Fluid Properties”

#### 4.2.8.3.1 Sub-category “Select Fluid”

Under this entry, you can select the thermal fluid used from a list.

## PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “Temperature Control”.
- Tap on the category “TAC/Manual”.
- Tap on the sub-category “Configuration auto”.
- Tap on the sub-category “Fluid Properties”.
- Tap on the sub-category “Select Fluid”.
- Select the thermal fluid used from the list.
- Confirm your choice by tapping on “OK”.

#### 4.2.8.3.2 Sub-category: “Bath/Circulation Volume”

Under this entry, you can enter the filling volume of the thermal fluid in your bath/circuit.

## PROCEDURE

- Go to the “Categories Menu”.
- Tap on the category “Temperature Control”.
- Tap on the category “TAC/Manual”.
- Tap on the sub-category “Configuration auto”.
- Tap on the sub-category “Fluid Properties”.
- Tap on the sub-category “Bath/Circulation Volume”.
- Read the message and confirm it by tapping on “OK”.
- Enter the filling volume using the number keypad that appears.
- Confirm your entry by tapping on “OK”.

4.2.8.3.3 Sub-category: "VPC/Bypass"

Under this entry, you can specify whether you use a bypass or not.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Fluid Properties".
- Tap on the sub-category "VPC/Bypass".
- Choose between the dialog entries "Bypass used" and "Bypass not used".
- Confirm your choice by tapping on "OK".

4.2.8.3.4 Sub-category: "Show Fluid"

This entry provides you with an overview of the settings entered.

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Fluid Properties".
- Tap on the sub-category "Show Fluid".
- Tap on "OK" after you have read/checked the entries.

4.2.8.4 Sub-category: "Display parameters"

Here you can display the set parameters in "Automatic mode".

PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration auto".
- Tap on the sub-category "Display parameters".
- Tap on "OK" after you have read/checked the parameters.

4.2.9 Sub-category: "Configuration manual"

**NOTE**

**Use of the "Expert mode" without a thorough knowledge of I&C technology.  
MATERIAL DAMAGE TO THE APPLICATION**

- Only use this mode if you have a thorough knowledge of I&C technology.

**INFORMATION**

In "Expert mode", "Configuration auto" is deactivated only "Configuration manual" is possible.

To enter settings in "Expert mode" requires in-depth knowledge of process measuring and control technology. Incorrect or insufficient settings can severely impair the function of the temperature control unit.

4.2.9.1 Sub-category: "Change parameters"

In this menu you manually configure the control parameters. If only the internal temperature is used as a target value, the control parameters are entered only under "Internal". If the process temperature is used as the target value, also the internal controller can be used, for example, when reaching the set point limit or in case of a Delta T limit. Consequently, parameter sets must be entered under all three items ("Internal", "Jacket" and "Process") when using the process temperature as a target value.

#### 4.2.9.1.1 Sub-category: "Internal"

Enter the new values for "KP", "Tn" and "Tv" here one after the other.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Change parameters".
- Tap on the sub-category "Internal".
- Enter the new "KP" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tn" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tv" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".

#### 4.2.9.1.2 Sub-category: "Jacket"

Enter the new value for "KP" here.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Change parameters".
- Tap on the sub-category "Jacket".
- Enter the new "KP" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Read the message and confirm it by tapping on "OK".

#### 4.2.9.1.3 Sub-category: "Process"

Enter the new values for "KP", "Tn" and "Tv" here one after the other.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Change parameters".
- Tap on the sub-category "Process".
- Enter the new "KP" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tn" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".
- Enter the new "Tv" value using the number keypad that appears.
- Confirm your entry by tapping on "OK".

#### 4.2.9.2 Sub-category: "Display parameters"

In this function, the set manual parameters are displayed.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".

- Tap on the sub-category "Configuration manual".
- Tap on the sub-category "Display parameters".
- Tap on "OK" after you have read/checked the parameters.

#### 4.2.9.3 Sub-category: "Control structure"

With this function, you have two different control structures available.

"Huber PID controller": Default setting

"Classic PID controller": This setting is exclusively used by Huber service engineers for service purposes.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Configuration manual".
- Tap on the sub-category "control structure".
- Choose between the dialogue entries "Huber PID controller" and "Classic PID controller".
- Confirm your choice by tapping on "OK".

#### 4.2.10 Sub-category: "Reset parameters"

With this function, you can reset the control parameters to the factory setting.

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Reset parameters".
- Read the message and confirm it by tapping on "Yes". The control parameters are reset/deleted. **The temperature control unit can only be operated again after a complete restart.**
- To do so, switch the temperature control unit off and back on again. The parameters have been reset.

#### 4.2.11 Sub-category: "Display parameters"

In this function, the set parameters are displayed. Depending on the previous setting, these will be the "Automatic control parameters" or the "Manual control parameters".

### PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "TAC/Manual".
- Tap on the sub-category "Display parameters".
- Tap on "OK" after you have read/checked the parameters.

#### 4.2.12 Freeze protection for temperature control unit

##### NOTE

**Non-compliance with the compatibility between the thermal fluid and your temperature control unit**

##### MATERIAL DAMAGE

- Observe the classification of your temperature control unit according to DIN 12876.
- Ensure the following materials are resistant with respect to the thermal fluid: Stainless steel 1.4301/ 1.4401 (V2A), copper, nickel, FKM, red bronze/brass, silver solder and plastic.
- The maximum viscosity of the thermal fluid must not exceed 50 mm<sup>2</sup>/s at the lowest working temperature!
- The maximum density of the thermal fluid may not exceed 1 kg/dm<sup>3</sup>!

**NOTE****Operation without freeze protection activated****COMPLETE DESTRUCTION OF THE TEMPERATURE CONTROL UNIT**

- When using water as a thermal fluid, freeze protection must be activated in the category “protective devices”.


Temperature control units with freeze protection option can use water as the only thermal fluid. Possible danger: Freezing of the heat exchanger at  $< 5\text{ °C}$ . An ice crystal (green flashing) indicates the activation of the freeze protection on the touch display. To find out whether your temperature control unit is equipped with a freeze protection, refer to the datasheet from page 73 in section »Annex«.

**PROCEDURE**

- Go to the “Categories Menu”.
- Tap on the category “Protection Options”.
- Tap on the category “Freeze Protection”.
- Either tap on the dialog entry “off” or “on” to turn the freeze protection on or off.
- Confirm your choice by tapping on “OK”.

**4.2.13 Setting the setpoint thresholds**

Overview of the temperature thresholds

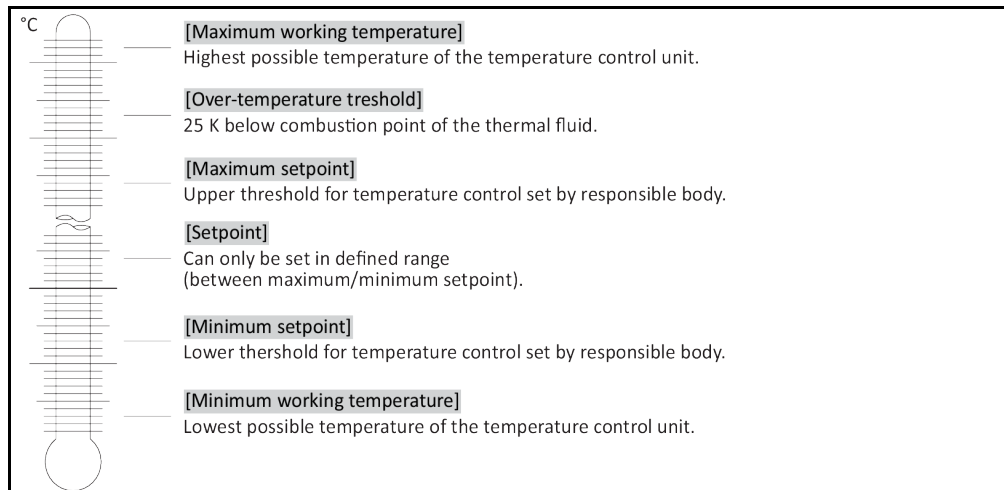

The limits for the minimum and maximum setpoints serve for the safety of your system. They should be set for the application field of the thermal fluid before starting the first temperature control and when changing the thermal fluid. The maximum setpoint limit limits the setpoint setting for the bath or flow temperature. The minimum setpoint limit protects against high viscosity or freezing at low temperatures. Setpoint settings are only available in the temperature range between the minimum and maximum setpoint limit.

**PROCEDURE**

- Go to the “Categories Menu”.
- Tap on the category “Protection Options”.
- Tap on the category “Setpoint Limits”.
- Tap on the sub-category “Minimum Setpoint”.
- Enter the new value, using the numeric keypad that appears.
- Confirm your entry by tapping on “OK”.
- In the display that follows, confirm your entry again by tapping on “OK”.  
The correct selection will be displayed graphically and the “Min. setpoint” will be changed promptly. If tapping on “OK” is not correct, it will be displayed graphically for 2 seconds. The display will then return to the category “Setpoint Limits”. Try changing the “Min. Setpoint” again.
- Tap on the sub-category “Maximum Setpoint”.
- Enter the new value, using the numeric keypad that appears.
- Confirm your entry by tapping on “OK”.
- In the display that follows, confirm your entry again by tapping on “OK”.

The correct selection will be displayed graphically and the “Max. Setpoint” will be changed promptly. If tapping on “OK” is not correct, it will be displayed graphically for 2 seconds. The display will then return to the category “Setpoint Limits”. Try changing the “Max. Setpoint” again.

**INFORMATION**

Check the set values for the minimum and maximum setpoint at any system change, especially when changing the thermal fluid.

4.2.14 Setting the setpoint

**PROCEDURE**

- Go to the “Home” screen.
- Tap on the keypad symbol next to “T<sub>Setpoint</sub>”.
- Enter a new setpoint using the number keypad that appears.  
**The following must apply:**  
 [Minimum setpoint limit] ≤ [Setpoint] ≤ [Maximum setpoint limit].  
 If these conditions are not met, a message will appear on the >Touchscreen< [88] and the entry will be ignored. In this case, delete the value just entered either with the “Arrow” key or with the “clear” key. Enter the setpoint again.
- Confirm your entry by tapping on “OK”.
- In the display that follows, confirm your entry again by tapping on “OK”.  
 The correct selection will be displayed graphically and the setpoint will be changed immediately. If tapping on “OK” is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the “Home” screen again. Try changing the setpoint again.

**4.3 Filling, venting, degassing and draining**

The illustration “connection diagram” can be found on page 73 in section »Annex«.

**CAUTION**

**Extremely hot / cold surfaces, connections and thermal fluids**

**BURNS/FREEZING OF LIMBS**

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

**NOTE**

**During an active circulation, the thermal fluid circuit is shut off by shut-off valves**

**MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT**

- Do not close the thermal fluid circuit during an active circulation by means of shut-off valves.
- Warm the thermal fluid to room temperature before stopping the circulation.

4.3.1 Filling, venting, degassing and draining the bath thermostat

**CAUTION**

**Non-compliance with the safety data sheet for the thermal fluid to be used**

**INJURIES**

- Risk of injury to the eyes, skin, respiratory tract.
- The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected.
- Observe the local regulations/work instructions.
- Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
- Danger of slipping because floor and work area are contaminated. Clean the work station and follow the instructions for the disposal of thermal fluid and material on page 15 in section »Proper disposal of resources and consumables«.

#### 4.3.1.1 Filling and venting the bath thermostat and the externally closed application

### PROCEDURE

- Check that a hose was installed at the **>Overflow< [12]** (if present). The other end of the hose must be inserted in a suitable collecting container. Excess thermal fluid will leak at this point when the temperature control unit is overfilled. Hose and container must be compatible with the thermal fluid and its temperature.
- Lift the **>Bath cover< [93]** from the temperature control unit.
- Carefully pour in suitable thermal fluid (see page 29 in section **»Information on the thermal fluids«**) using the filling accessories (funnel and/or beaker). While filling, pay attention to any necessary measures such as grounding the containers, funnels and other accessories. The thermal fluid can flow via the hose connection to the external application. Follow the instructions for cleaning the filling accessories on page 15 in section **»Proper disposal of resources and consumables«**.
- Switch on the temperature control unit.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Start venting".
- Confirm your choice by tapping on "OK".
- Using the number keypad that appears, enter the time interval for venting. The default is 0.5 minutes.
- Confirm your entry by tapping on "OK" to start the filling process. The filling process is complete when the bath vessel is filled sufficiently. For temperature control units with heater, each heating coil of the heater must be below the thermal fluid level. This corresponds to the minimum level in the bath/minimal filling. There must be no air bubbles left in the external application (e.g. in the double jacket). The control electronics monitors the fill quantity and displays the level on the graphic display.
- Run the "venting program" for a few minutes. This ensures that trapped air bubbles can escape. These bubbles would result in a safety shutdown during an actual temperature control process.
- Stop venting. To do this, go to the category "Temperature Control".
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Stop venting".
- Confirm your choice by tapping on "OK". Venting is stopped and the pump continues to run for approx. 30 seconds. Wait until the pump stops.
- Put the **>bath cover< [93]** back onto the bath opening.

#### INFORMATION

The volume expansion of the thermal fluid depends on the working temperature range you wish to work in. Do not go below the minimum bath level/minimal level when working at the "lowest" working temperature and there should be no overflow from the expansion vessel/temperature control unit when working at the "highest" working temperature. In case of overfilling, drain the excess amount of thermal fluid (see page 56 in Section **»Draining the bath thermostat«**). Check if the thermal fluid can be reused. On page 15 observe Section **»Proper disposal of resources and consumables«**.

- Check the level in the collecting container. Empty the container when necessary and dispose of its contents properly.
- In case of overfilling, drain thermal fluid via the **>Drain< [8]** into a suitable container (see page 56 the section **»Draining the bath thermostat«**). Check if the thermal fluid can be reused. Please note on page 15 the section **»Proper disposal of resources and consumables«**.

#### INFORMATION

**Venting and degassing** must be performed especially during commissioning and after changing the thermal fluid. This is the only way to ensure trouble-free operation. After venting is completed, follow the steps on page 55 in section **»Degassing of bath thermostat«**

#### 4.3.1.2 Degassing of bath thermostat

#### CAUTION

##### Hot or cold thermal fluid and surfaces BURNS TO LIMBS

- Avoid direct contact with the thermal fluids or the surfaces.
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).

**INFORMATION**

When changing from low-boiling thermal fluid (low-boiling components) to higher boiling thermal fluids, remains of the low-boiling component may remain in the temperature control unit. Depending on the working temperature, the low-boiling component begins to boil and gas bubbles are formed that cause the pump pressure to momentarily collapse. This may cause a safety shutdown. The gas bubbles reach the bath opening and can escape.

If ice crystals form on the evaporator coil, water has accumulated in the thermal fluid. Degas if this is the case, to avoid damage to the temperature control unit.

Thermal fluids are more or less hygroscopic (water-attracting). This effect increases, the lower the working temperature. The de-gassing mode below, which must be **permanently monitored**, also helps you remove any water residues from the temperature control circuit.

**PROCEDURE**

- Follow venting with the degassing operation. Prerequisite: You have filled the temperature control unit in accordance with the instructions on page 55 in Section »Filling and venting the bath thermostat and the externally closed application« and/or cleaned it as per page 66 in Section »Rinsing the thermal fluid circuit«.
- Go to the "Home" screen.
- Tap on the keypad symbol next to "T<sub>Setpoint</sub>".
- Enter a setpoint using the number keypad that appears. This setpoint must be below the lower boiling thermal fluid. This setpoint will be increased in 10 K steps during the degassing process up to the maximum working temperature.
- Confirm your entry by tapping on "OK".
- In the display that follows, confirm your entry again by tapping on "OK". The correct selection will be displayed graphically and the setpoint will be changed immediately. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will return to the "Home" screen. Try changing the setpoint again.
- Tap on the "Start" touch button.
- Confirm the start of temperature control by tapping on "OK". The correct selection will be displayed graphically and temperature control will start immediately. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will return to the "Home" screen. Try starting the thermoregulation again.
- Carry out temperature control to the entered setpoint until no more gas bubbles rise up.
- Increase the setpoint by 10 K and carry out temperature control until no more gas bubbles rise up.
- Repeat increasing the setpoint by 10 K until the maximum working temperature of the thermal fluid has been reached.
- Tap on the "Stop" touch button as soon as no more gas bubbles rise up at the thermal fluid's maximum working temperature.
- Confirm the stop of the temperature control process by tapping on "OK". The correct selection will be displayed graphically. Thermoregulation is stopped immediately and the pump continues to run for approx. 30 seconds. Wait until the pump stops. If tapping on "OK" is not correct, it will be displayed graphically for 2 seconds. The display will return to the "Home" screen. Try stopping the temperature control process again. The de-gassing process is complete.

**4.3.1.3 Draining the bath thermostat****CAUTION****Hot or very cold thermal fluid****SEVERE BURNS/FROSTBITE OF LIMBS**

- Before draining, ensure that the thermal fluid has room temperature (20 °C).
- If, at this temperature, the thermal fluid is too viscous to be drained: Control the temperature of the thermal fluid for a few minutes until the viscosity will allow drainage. Never control the temperature of the thermal fluid when the Drain is open.
- Danger of burns when draining thermal fluids at temperatures above 20 °C.
- Wear your Personal Protective Equipment for draining.
- Only drain using suitable drainage hose and collecting container. These must be compatible with the thermal fluid and its temperature.



## PROCEDURE

### Baths with >Drain valve< [4]

- Remove the knurled screw at the >Drain< [8].
- Connect a suitable drain hose to the >Drain< [8].
- Place the other end of the hose in a suitable container.
- Open the >Drain valve< [4] by turning it counterclockwise (turn 90° left as far as it will go). The thermofluid will flow from the external application via the bath vessel and the draining hose into the container. Check if the thermofluid can be reused. On page 15 observe section »**Proper disposal of resources and consumables**«.
- Wait until the external application and the bath are empty.
- Open the connection >Circulation flow< [1].
- Open the connection >Circulation return< [2].
- Leave the temperature control unit open for a while to allow it to dry out and the residue to drain. Without screw caps and with the >Drain valve< [4] open.
- Close the >Drain valve< [4] by turning it clockwise (turn 90° right as far as it will go).
- Close the connection >Circulation flow< [1].
- Close the connection >Circulation return< [2].
- After drying out, remove the drain hose and re-fit the knurled screw to the >Drain< [8].
- The bath is now drained.

### Baths without >Drain valve< [4]

- Have a suitable container ready to catch the thermal fluid.
- Open the knurled screw at the >Drain< [8]. As soon as you have opened the knurled screw, the thermal fluid will flow from the external application over the bath and into the container. Wait until the external application and the bath are empty. Check if the thermal fluid can be reused. On page 15 observe Section »**Proper disposal of resources and consumables**«.
- Open the connection >Circulation flow< [1].
- Open the connection >Circulation return< [2].
- Leave the temperature control unit open for a while for the residue to fully drain and to allow it to dry out (without screw caps).
- Close the connection >Circulation flow< [1].
- Close the connection >Circulation return< [2].
- Re-fit the knurled screw to the >Drain< [8].
- The bath is now drained.

## 5 Normal operation

### 5.1 Automatic operation

#### CAUTION

**Extremely hot / cold surfaces, connections and thermal fluids**

**BURNS/FREEZING OF LIMBS**

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

#### NOTE

**During an active circulation, the thermal fluid circuit is shut off by shut-off valves**

**MATERIAL DAMAGE TO THE CIRCULATING PUMPS INSTALLED IN THE TEMPERATURE CONTROL UNIT**

- Do not close the thermal fluid circuit during an active circulation by means of shut-off valves.
- Warm the thermal fluid to room temperature before stopping the circulation.

#### 5.1.1 Temperature control

##### 5.1.1.1 Starting the temperature control process

The temperature control process can be started after filling and complete venting.

### PROCEDURE

- Go to the “Home” screen.
- Tap on the “Start” touchbutton.
- Confirm the start of temperature control by tapping on “OK”.  
The correct selection will be displayed graphically and temperature control will start immediately. If tapping on “OK” is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the “Home” screen again. Try to start the temperature control unit again.

##### 5.1.1.2 Ending the temperature control process

#### NOTE

**When the temperature control unit is switched off, the thermal fluid temperature is higher/lower than room temperature**

**DAMAGE TO THE TEMPERATURE CONTROL UNIT AND THE GLASS APPARATUS/APPLICATION**

- Bring the thermal fluid up to room temperature using the temperature control unit.
- Do not close the shut-off valves in the thermal fluid circuit.

The temperature control can be ended at any time, the pump continues to run for approx. 30 seconds. The compressor is switched off after the stepper motor valve has reached a defined position to control the cold output.

### PROCEDURE

- Go to the “Home” screen.
- Tap on the “Stop” touchbutton.
- Confirm the stop of temperature control by tapping on “OK”.  
The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on “OK” is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the “Home” screen again. Try to stop the temperature control unit again.

#### INFORMATION

The compressor is not switched off until the stepper motor valve has reached a defined position. The status line [Field 10] displays the relevant information.

## 5.1.2 Temperature control via a created temperature control program

### 5.1.2.1 Starting the temperature control program

A temperature control program can be started after filling and complete venting.

## PROCEDURE

- Go to the "Categories Menu".
- Tap on the category "Programmer/Ramp".
- Tap on the category "Start/stop Program".
- Tap on the dialog entry of the temperature control program to be started.
- Confirm your choice by tapping on "OK".
- Read the message and confirm it. You temperature control unit will start the temperature control program and the temperature control programmed in it will start.
- Read the Note and confirm by tapping on "OK".

### 5.1.2.2 Ending/cancelling the temperature control program

#### NOTE

**When the temperature control unit is switched off, the thermal fluid temperature is higher/lower than room temperature**

#### **DAMAGE TO THE TEMPERATURE CONTROL UNIT AND THE GLASS APPARATUS/APPLICATION**

- Bring the thermal fluid up to room temperature using the temperature control unit.
- Do not close the shut-off valves in the thermal fluid circuit.

You can either end the temperature control automatically using the parameters predefined in the temperature control program or end/interrupt the temperature control at any time manually. The temperature control is switched off immediately after, the pump continues to run for approx. 30 seconds. The compressor is switched off after the stepper motor valve has reached a defined position to control the cold output.

#### Manual ending/cancelling

## PROCEDURE

- Go to the "Home" screen.
- Tap on the "Stop" touchbutton.
- Confirm the stop of temperature control by tapping on "OK".  
The correct selection will be displayed graphically and temperature control will stop immediately and the pump with overrun for about 30 seconds. Wait until the pump stops. If tapping on "OK" is not correct, this is displayed graphically for 2 seconds. After this, the display will return to the "Home" screen again. Try to stop the temperature control unit again.

#### INFORMATION

The compressor is not switched off until the stepper motor valve has reached a defined position. The status line [Field 10] displays the relevant information.

## 6 Interfaces and software update

**NOTE**

**Connecting to the interfaces at the temperature control unit during operation**

**DAMAGE TO THE INTERFACES**

- When devices in operation are connected with interfaces of the temperature control unit, interfaces may get damaged.
- Before connecting, ensure the temperature control unit and the device to be connected are turned off.

**NOTE**

**The specifications of the interface used are not being met.**

**PROPERTY DAMAGE**

- Only connect components that meet the specifications of the interface used.

**INFORMATION**

The use of PB commands is described in our “Data communications PB” manual. This manual can be downloaded from [www.huber-online.com](http://www.huber-online.com).

### 6.1 Interfaces at the “Pilot ONE®” controller


**NOTE**

**The Pilot ONE controller is not operated behind a firewall**

**PROPERTY DAMAGE**

- Operate the controller Pilot ONE exclusively behind a firewall, if the local subnet is connected to the Internet or another risk-prone network.
- The best available technology is to be applied in order to provide sufficient security for the LAN!

Standard interfaces at the “Pilot ONE” top side

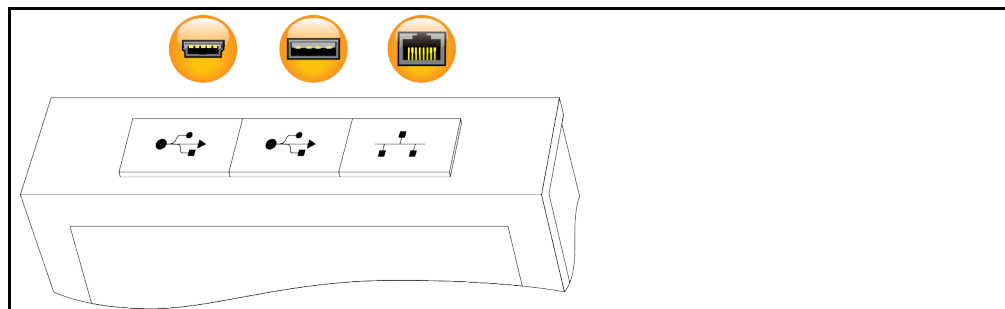

#### 6.1.1 10/100 Mbps Ethernet for RJ45 network sockets



This is a fast and flexible interface. Standard 10/100 Mbps interface (Fast Ethernet), can be connected to any existing Ethernet network. Because this interface can also be connected to very large networks, the IT “Best Practices” (firewall) must be observed.

**Usage:**

Also - to be able to communicate with the “Pilot ONE” controller - the communication enable must be issued. This is an additional safety feature that prevents persons - possibly unintentionally - connecting to the wrong machine and implementing incorrect temperature control specifications. The following restrictions are possible:

- Deactivated
- Always on (PLC)
- 12h Inactivity Timer
- 10min Inactivity Timer

If, for example, “10min Inactivity Timer” is selected, the connection must be made within 10 minutes after confirmation at the control. If this does not happen, the connection is refused.

**INFORMATION** Communication with the Pilot ONE is via TCP (Transmission Control Protocol), Port 8101. The interfaces used must meet the specifications of the generally accepted standards.

### 6.1.2 USB-2.0 interface

**INFORMATION** The interfaces used must meet the specifications of the generally accepted standards.

#### 6.1.2.1 USB-2.0 interface, host



USB-2.0 connection (for connector A), e.g. for data memories.

#### 6.1.2.2 USB-2.0 interface, device




USB-2.0 connection (for Mini-B connector) for communicating with a computer.

## 6.2 Interfaces on the temperature control unit

### 6.2.1 Interfaces at the back

Interfaces at the back of the immersion circulator (exemplary layout)

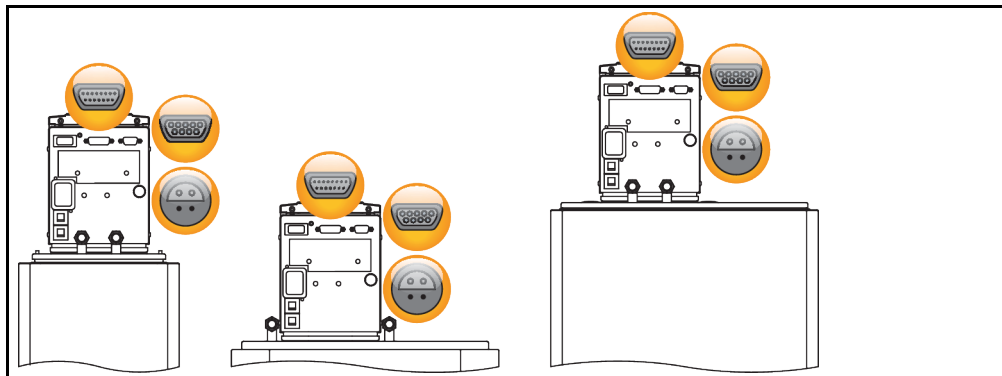

#### 6.2.1.1 Service interface



This interface is exclusively used by Huber service engineers for service purposes. An adapter cable makes this interface a RS232 serial port.

#### 6.2.1.2 Female RS232 serial



A PC, a SPS or a Process Control System (PCS) can be connected to this jack for remote control of the controller electronics. Before plugging in the cable, check the settings in the “Interfaces” category and adjust if necessary.

**INFORMATION** The interfaces used must meet the specifications of the generally accepted standards.

Pin assignment (front view)



Pin assignment

Pin	Signal	Description
2	RxD	Receive Data
3	TxD	Transmit Data
5	GND	Signal GND

**6.2.1.3 Connection jack for Pt100 process controller sensor**



A temperature sensor located in the connected application (Pt100, 4-wire technology, Lemosa connector) is connected to the Pt100 connection jack. The external actual temperature is then recorded and the operating temperature of the temperature control unit is permanently calculated and adjusted.

**INFORMATION**

Depending on the operating temperature, isolation losses and exothermic heat, the operating temperature (flow temperature) of the application can be significantly less than the setpoint of the application. In this context, safety-critical thresholds for the temperature control fluid must be strictly observed.

The control results contained in the data sheet can only be achieved with **shielded** sensor leads. We recommend the external Pt100 process control sensor from the Huber accessories program.

Pin assignment (front view)



Pin assignment

Pin	Signal
1	I+
2	U+
3	U-
4	I-

Pt100

Pin 1: I+    Pin 2: U+    Pin 3: U-    Pin 4: I-

**6.3 Firmware update**


An instruction for running a firmware update can be found at [www.huber-online.com](http://www.huber-online.com).

## 7 Service/maintenance

### 7.1 Messages from the temperature control unit

Messages output by the temperature control unit can be divided into various classes.

Follow the instructions displayed on the >Touchscreen< [88]. Once a message has been acknowledged, a symbol is output on the >Touchscreen< [88]. Tapping the symbol takes you to an overview of all messages in chronological order.

Displayed symbol: 

### 7.2 Electrical fuse

The thermal overcurrent circuit breakers for all-pole disconnection (L and N) are located at the back. In case of a fault (no function and/or no display) please first check if the overcurrent circuit breakers have tripped. If the overcurrent circuit breakers trigger again immediately after reversing, please unplug the power cord and contact Customer Support immediately (the phone number can be found on page 72 in section »Contact data«).

### 7.3 Replacing the “Pilot ONE®”




**Replacing the electronics while the temperature control unit is operating**

**MORTAL DANGER FROM FIRE**

- Stop an ongoing temperature control process.
- Turn off the temperature control unit.
- Also disconnect the temperature control unit from the power supply.

In case of fault, you can replace the “Pilot ONE” yourself. If you have questions or difficulties, please contact your dealer, your local office or our Customer Support.

Replacing the  
“Pilot ONE”

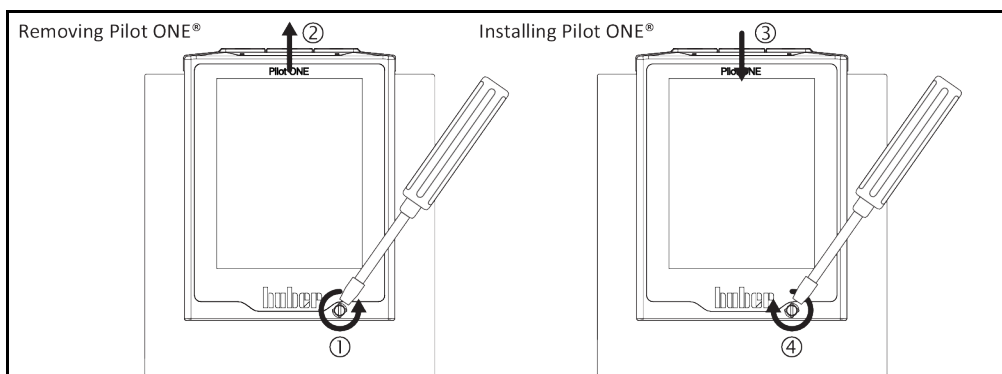

## PROCEDURE

- Turn off the temperature control unit.
- Disconnect the temperature control unit from the power supply.
- Release the >Fixing screw for Pilot ONE< [89] at the front of the housing.
- Carefully pull the “Pilot ONE” away upwards.
- Carefully insert the replacement “Pilot ONE”.
- Close the >Fixing screw for Pilot ONE< [89] at the front of the housing.
- Connect the temperature control unit to the power supply.
- Turn on the temperature control unit.

## 7.4 Maintenance



### Cleaning/maintenance while the temperature control unit is operating

#### MORTAL DANGER FROM ELECTRIC SHOCK

- Stop an ongoing temperature control process.
- Turn off the temperature control unit.
- Also disconnect the temperature control unit from the power supply.

### NOTE

### Performing maintenance work not described in these operation manual

#### MATERIAL DAMAGE ON THE TEMPERATURE CONTROL UNIT

- Please contact Huber for maintenance work that is not described in these operation manual.
- Maintenance work not described in these operation manual is reserved for qualified specialists trained by Huber.
- Safety-relevant components may only be replaced by equivalent ones. The specified safety values for the respective component must be observed.

### 7.4.1 Function check and visual inspection

Control intervals

Cooling*	Description	Maintenance interval	Comment	Person responsible
A/W	Visually inspect hoses and hose connections	Prior to switching on the temperature control unit	Exchange leaking hoses and hose connections prior to switching on the temperature control unit. Please see on page 65 the section »Replacing temperature control or coolant hoses«.	Operating company and / or operators
A/W	Check the fill level in the collecting container at the <b>&gt;Overflow&lt;</b> [12] (if present)	Prior to switching on the temperature control unit	Check the fill level in the collecting container. Empty as required. Follow the instructions for the disposal of thermal fluid on page 15 in the section »Proper disposal of resources and consumables«.	Operating company and / or operators
A/W	Check the power supply cable	Prior to switching on the temperature control unit or on relocation	Do not start the temperature control unit if the power supply cable is damaged.	Qualified electrician (BGV A3)
A/W	Thermal fluid inspection	As required	–	Operating company and / or operators
A/W	Overtemperature protection (OT) - functional check	Every month or after changing the thermal fluid.	Please see on page 42 the section »Testing overtemperature protection for functionality«	Operating company and / or operators
A/W	Check the temperature control unit for damage and stability	Every 12 months or after a change of location	–	Operating company and / or operators
W	Check the cooling water quality	Every 12 months	Descale the cooling water circuit as required. Documentation on water quality is available at: <a href="http://www.huber-online.com">www.huber-online.com</a>	Operating company and / or operators
A/W	Exchange safety-relevant electric and electromechanical components	20 years	Have the exchange only carried out by certified personnel (such as Huber service engineers). Please contact Customer Support. The telephone number can be found on page 72 in the section »Contact data«	Operating company

\*A = Air cooling; W = Water cooling; U = Applicable only for Unistats



## 7.4.2 Replacing temperature control or coolant hoses

Replace defective temperature control and/or coolant hoses **before** turning on the temperature control unit.

### 7.4.2.1 Replacing temperature control hoses

#### PROCEDURE

- Drain the temperature control unit as described on page 56 in Section »**Draining the bath thermostat**«.
- Replace defective temperature control hoses. When disposing of them, observe on page 15 Section »**Proper disposal of resources and consumables**«.
- Reconnect your external application as described on page 26 in Section »**Connecting externally closed application**«.
- Fill the temperature control unit with thermal fluid as described on page 55 in Section »**Filling and venting the bath thermostat and the externally closed application**«.
- Vent the temperature control unit as described on page 55 in Section »**Filling and venting the bath thermostat and the externally closed application**«.
- Restart the temperature control unit in normal mode.

### 7.4.2.2 Replacing coolant hoses

#### PROCEDURE

- Drain the cooling water as described on page 70 in section »**Draining the cooling water**«.
- Replace the defective coolant hoses. When disposing of them, observe on page 15 section »**Proper disposal of resources and consumables**«.
- Reconnect the temperature control unit to the building's cooling water supply. Proceed as described on page 24 in section »**Temperature control units with cooler**«.
- Restart the temperature control unit in normal mode.

## 7.5 Thermal fluid inspection, replacement and circuit cleaning

The illustration "connection diagram" can be found on page 73 in section »**Annex**«.

### CAUTION

#### Extremely hot / cold surfaces, connections and thermal fluids

##### BURNS/FREEZING OF LIMBS

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

### 7.5.1 Thermal fluid inspection

### CAUTION

#### Thermal fluid is not inspected on a regular basis

##### BURNS DUE TO REDUCED BOILING POINT

- Regularly check your thermal fluid whether it meets the specifications in the safety data sheet.

### NOTE

#### Thermal fluid is not inspected on a regular basis

##### DAMAGE TO THE HEAT EXCHANGER AND/OR ELECTROMECHANICAL PARTS.

- Regularly check your thermal fluid whether it meets the specifications in the safety data sheet.

**INFORMATION**

**Oxidation**

Oxidation ages the thermal fluid and change its characteristics (e.g. a reduced boiling point). When controlling high temperatures, a reduced boiling point may cause overflow of very hot thermal fluids. It may cause serious burns of the limbs.

**Hygroscopy**

When continuously thermoregulating below room temperature, hygroscopy causes the thermal fluid to accumulate water in the course of time. Such a liquid mixture causes the evaporator to burst when thermoregulating in the minus range. This is caused by the water in the liquid mixture, which forms ice crystals on the evaporator. When thermoregulating high temperatures with such a liquid mixture, the boiling point is reduced. When controlling high temperatures, a reduced boiling point may cause overflow of very hot thermal fluids. It may cause serious burns of the limbs. Hygroscopy can change the mixing ratio of a water-ethylene-glycol mixture.

7.5.2 Rinsing the thermal fluid circuit

**DANGER**

**Setpoint and overtemperature protection are not adjusted to the thermofluid**

**MORTAL DANGER FROM FIRE**

- The cut-out value of the overtemperature protection **must** be adapted to the thermofluid. Set the cut-out value of the overtemperature protection 25 K below the fire point of the thermofluid.
- The setpoint set during rinsing **must** be adjusted to the thermofluid used.

**CAUTION**

**Non-compliance with the safety data sheet for the thermal fluid to be used**

**INJURIES**

- Risk of injury to the eyes, skin, respiratory tract.
- The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected.
- Observe the local regulations/work instructions.
- Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
- Danger of slipping because floor and work area are contaminated. Clean the work station and follow the instructions for the disposal of thermal fluid and material on page 15 in section »Proper disposal of resources and consumables«.


**NOTE**

**Mixing different thermofluids in a thermal fluid circuit**

**PROPERTY DAMAGE**

- Do **not** mix different types of thermofluid (such as mineral oil, silicone oil, synthetic oil, water, etc.) in a thermofluid circuit.
- The thermal fluid circuit **must** be rinsed when changing from one type of thermal fluid to another. No residues of the previous type of thermal fluid may remain in the thermal fluid circuit.

Example: Connecting a short circuit hose

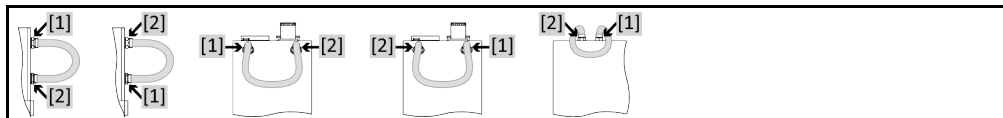

The inner components of the temperature control unit must be dried out. Need to avoid boiling retardation during future uses (e.g. use of a silicone oil at temperatures above about 100 °C).

**PROCEDURE**

- Drain the temperature control unit as described on page 56 in section »Draining the bath thermostat«.

**INFORMATION**

Residual thermal fluid can remain in the pump chamber and the internal lines after draining. Leave the temperature control unit with open valves for a while.

- Leave the drain hose mounted to the >Drain< [8].
- Check the fill level in the collecting container at the end of the drain hose. Follow the instructions on page 15, section »Proper disposal of resources and consumables« for the disposal of thermofluid.
- Close the drain valves on the temperature control unit by turning them clockwise (turn 90° right as far as it will go).

- Connect the >Circulation flow< [1] with the >Circulation return< [2] to the temperature control using a bypass hose.

**INFORMATION**

Perform the following steps without attaching a short circuit hose, if the application used by you (externally closed) is also dirty. In this case, leave your externally closed application connected to the temperature control unit. This rinses the temperature control unit and your application at the same time.

- **Fill** the system (minimum fill level) with the thermal fluid you wish to use. The filling process is described on page 55 in section »Filling and venting the bath thermostat and the externally closed application«.
- **Vent** the system as described on page 55 in section »Filling and venting the bath thermostat and the externally closed application«.
- Adjust the **setpoint** and the cut-out value of the **overtemperature protection** to the thermal fluid used. The procedure can be found on page 54 in section »Setting the setpoint« and on page 41 in section »Setting the overtemperature protection«.
- Go to the "Categories Menu".
- Tap on the category "Temperature Control".
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Start temperature control process".
- Confirm your choice by tapping on "OK". The length of rinsing depends on the level of soiling.
- Tap on the category "Start/Stop".
- Tap on the dialog entry "Stop temperature control process".
- Confirm your choice by tapping on "OK". The temperature control process is stopped.
- **Drain** the temperature control unit as described on page 55 in section »Filling and venting the bath thermostat and the externally closed application«.
- Repeat the steps "Filling", "Venting", "Start/Stop temperature control process" and "Draining" until the drained thermal fluid remains clear.
- Remove the bypass hose after completely emptying the temperature control unit.

**INFORMATION**

Leave an application connected, if you have simultaneously rinsed a used application (externally closed).

- Leave the >Drain< [8] open for a while to allow the thermofluid to evaporate in the temperature control unit.
- Close the >Drain< [8] once the residual thermofluid has evaporated.
- Dismount the drain hose.
- Remove the collecting container.
- Discard the collecting container, including its contents, properly. Follow the instructions on page 15, section »Proper disposal of resources and consumables« for the disposal of thermofluid.
- Re-connect your application. (Only if you have rinsed the thermofluid circuit using a bypass hose.)
- Fill the temperature control unit with thermofluid as described on page 55 in section »Filling and venting the bath thermostat and the externally closed application«.
- Drain the temperature control unit as described on page 55 in section »Filling and venting the bath thermostat and the externally closed application«. An externally open application does not need to be vented.
- Start the "degassing" function as described on page 55 in section »Degassing of bath thermostat«. An externally open application does not need to be de-gassed.
- Restart the temperature control unit in normal mode.

## 7.6 Cleaning the surfaces

 **CAUTION**

### Extremely hot / cold surfaces, connections and thermal fluids

#### BURNS/FREEZING OF LIMBS

- Surfaces, connections and tempered thermal fluids can be extremely hot or cold depending on the operating mode.
- Avoid direct contact with surfaces, connections and thermal fluids!
- Wear your personnel protective equipment (e.g. temperature-resistant safety gloves, safety goggles).

**NOTE**

**Exposed plug contacts**

**DAMAGE CAUSED BY FLUID INGRESS**

- Protect unused plug contacts with the protective caps supplied.
- Clean surfaces only with a damp cloth.

A standard stainless steel cleaning agent is suitable for cleaning the stainless steel surfaces. Carefully clean painted surfaces (damp only) using a solution of sensitive-fabrics detergent. Follow the instructions on page 15, section »**Proper disposal of resources and consumables**« for the disposal of cleaning agents and material.

## 7.7 Plug contacts

**NOTE**

**Exposed plug contacts**

**DAMAGE CAUSED BY FLUID INGRESS**

- Protect unused plug contacts with the protective caps supplied.
- Clean surfaces only with a damp cloth.

Protective caps are supplied for all plug contacts. Make sure that any plug contacts not required are protective with the caps.

## 7.8 Decontamination/repairs

 **CAUTION**

**Returning a not decontaminated temperature control unit for repair**

**PHYSICAL INJURY AND PROPERTY DAMAGE CAUSED BY HAZARDOUS MATERIALS IN OR ON THE TEMPERATURE CONTROL UNIT**

- Carry out appropriate decontamination.
- The decontamination process depends on the type and quantity of the materials used.
- Consult the relevant safety data sheet.
- You will find a prepared return receipt at [www.huber-online.com](http://www.huber-online.com).

As the responsible body you are responsible for carrying out decontamination **before** third-party personnel come into contact with the temperature control unit / accessory. Decontamination must be carried out **before** the temperature control unit / accessory is returned for repair or inspection. Attach a clearly visible written notice stating that the temperature control unit / accessory has been decontaminated.

To simplify the process, we have prepared a form for you. This is available for download at [www.huber-online.com](http://www.huber-online.com).

## 8 Shutting down

### 8.1 Safety instructions and basic principles

#### DANGER

**Connection/adjustment to the power supply not carried out by an electrician and/or connection to a power socket without protective earth (PE)**

##### **MORTAL DANGER FROM ELECTRIC SHOCK**

- Have the connection/adjustment to the power supply carried out by an electrician.
- Always connect the temperature control unit to safety sockets (PE).

#### DANGER

**Damaged power cable/power cable connection**

##### **MORTAL DANGER FROM ELECTRIC SHOCK**

- Do not start up the temperature control unit.
- Isolate the temperature control unit from the power supply.
- Have the power supply cable/power supply connection replaced and inspected by an electrician.
- Do not use a power cable that is longer than **3 m**.

#### WARNING

**Risk of tipping due to unstable temperature control unit**

##### **SERIOUS INJURY AND PROPERTY DAMAGE**

- Avoid risk of tipping due to unstable temperature control unit.

#### CAUTION

**Non-compliance with the safety data sheet for the thermal fluid to be used**

##### **INJURIES**

- Risk of injury to the eyes, skin, respiratory tract.
- The safety data sheet for the thermal fluid to be used must be read prior to using it and its content must be respected.
- Observe the local regulations/work instructions.
- Wear your personal protective equipment (e.g. temperature-resistant safety gloves, safety goggles, safety footwear).
- Danger of slipping because floor and work area are contaminated. Clean the work station and follow the instructions for the disposal of thermal fluid and material on page 15 in section »Proper disposal of resources and consumables«.

#### CAUTION

**Hot or very cold thermal fluid**

##### **SEVERE BURNS/FROSTBITE OF LIMBS**

- Before draining, ensure that the thermal fluid has room temperature (20 °C).
- If, at this temperature, the thermal fluid is too viscous to be drained: Control the temperature of the thermal fluid for a few minutes until the viscosity will allow drainage. Never control the temperature of the thermal fluid when the Drain is open.
- Danger of burns when draining thermal fluids at temperatures above 20 °C.
- Wear your Personal Protective Equipment for draining.
- Only drain using suitable drainage hose and collecting container. These must be compatible with the thermal fluid and its temperature.

#### **INFORMATION**

All safety instructions are important and must be followed accordingly during working operations!

### 8.2 Switch-off

#### PROCEDURE

- Turn off the temperature control unit.
- Disconnect the temperature control unit from the power supply connection.

### 8.3 Draining the temperature control unit

#### PROCEDURE

- Drain the temperature control unit as described from page 54 in section »Filling, venting, degassing and draining«.

### 8.4 Draining the cooling water

#### INFORMATION

This section must be observed when using water-cooled temperature control units.

#### 8.4.1 Draining process

#### CAUTION

##### Pressurized cooling water connections

##### RISK OF INJURY

- Wear your personnel protective equipment (e.g. safety goggles).
- Carefully open the cooling water connection. Open slowly (1-2 signal edges) and drain the cooling water slowly.

#### NOTE

##### Building side isolating valves are not closed

##### DAMAGE BY ROOM FLOODING

- Close the building's isolating valves in the cooling water supply and return lines.

#### PROCEDURE

- Close the building's isolating valve in the cooling water supply and return lines.
- Place the collecting container below the input and output of the >Cooling coil< [29].
- Unscrew the connection at the >Cooling coil< [29]. The cooling water will begin to drain from the lines.
- Remove the cooling water from the >Cooling coil< [29]. Allow the cooling water to fully drain to prevent the risk of freezing during transport and storage!

### 8.5 Deinstalling the collecting container

#### PROCEDURE

- Remove the hose from the collecting container.
- Follow the instructions for the disposal of thermal fluid on page 15 in Section »Proper disposal of resources and consumables«.
- Dismantle the hose at the >Overflow< [12].

### 8.6 Uninstalling an external application

#### PROCEDURE

- Disconnect the external application from the temperature control unit.

### 8.7 Packing

Use the original packaging wherever possible! Further information can be found on page 20 in section »Unpacking«.

## 8.8 Shipping


### NOTE

#### Temperature control unit transported incorrectly

##### PROPERTY DAMAGE

- Do not transport by truck on the castors or feet.
- Comply all requirements in this section to avoid damage to the temperature control unit.

Pallet with squared timber for free-standing units

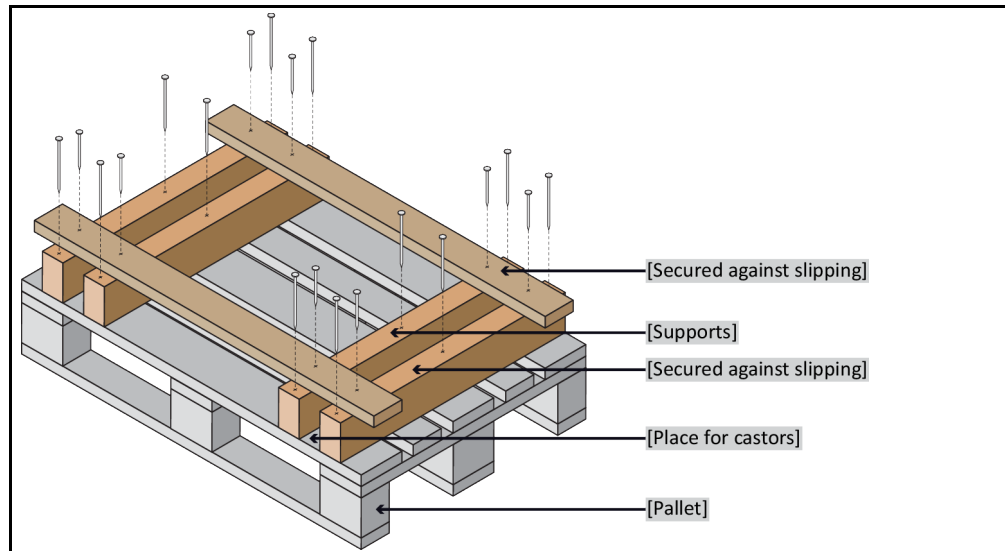

Transport using the lugs, if fitted, on the top of the temperature control unit. Do not transport the temperature control unit alone and without aids.

- Always use the original packaging for transportation.
- Indicate the upright transport position with arrows on the packaging.
- Always transport the temperature control unit upright on a pallet!
- Protect attachments from damage during transportation!
- During transport, place the temperature control unit on squared timber to protect the castors/feet.
- Secure with tensioning belts/lashing straps that are suitable for the weight.
- Additionally secure (depending on model) with plastic film, cardboard and straps.

## 8.9 Disposal

### NOTE

#### Improper disposal

##### ENVIRONMENTAL DAMAGE

- Spilled/leaked thermofluid must be discarded immediately and correctly. Follow the instructions for the disposal of thermofluid and material on page 15 in section »Proper disposal of resources and consumables«.
- To avoid environmental damage, have “disused” temperature control units disposed of exclusively by approved waste management companies (e.g. refrigeration and air conditioning companies).

Huber temperature control units and Huber accessories are made of high quality, recyclable materials. For example: Stainless steel 1.4301 / 1.4401 (V2A), copper, nickel, FKM, Perbunan, NBR, ceramic, carbon, Al-Oxid, red brass, brass, nickel-plated brass and silver solder. Proper recycling of the temperature control unit and accessories can actively help reduce CO<sub>2</sub> emissions in the production of these materials. Follow the laws and regulations of your jurisdiction when disposing material.

## 8.10 Contact data

### INFORMATION

Contact your supplier or local specialist retailer **prior** to returning the temperature control unit. The contact data can be found on our homepage [www.huber-online.com](http://www.huber-online.com) under the heading „Contact“. Please keep the serial number of the temperature control unit ready. The serial number can be found on the nameplate of the temperature control unit.

### 8.10.1 Telephone number: Customer Support

If your country is not mentioned in the list below: The responsible service partner can be found on our homepage [www.huber-online.com](http://www.huber-online.com) under the heading „Contact“.

- Huber Deutschland: +49 781 9603 244
- Huber China: +86 (20) 89001381
- Huber India: +91 80 2364 7966
- Huber Ireland: +44 1773 82 3369
- Huber Italia: +39 0331 181493
- Huber Swiss: +41 (0) 41 854 10 10
- Huber UK: +44 1773 82 3369
- Huber USA: +1 800 726 4877 | +1 919 674 4266

### 8.10.2 Telephone number: Sales

Telephone: +49-781-9603-123

### 8.10.3 Email address: Customer Support

Email: [support@huber-online.com](mailto:support@huber-online.com)

## 8.11 Certificate of Compliance

Please read page 68, section »Decontamination/repairs«.



## 9 Annex

# Inspired by **temperature** designed for you

Peter Huber Kältemaschinenbau AG  
Werner-von-Siemens-Str. 1  
77656 Offenburg / Germany

Telefon +49 (0)781 9603-0  
Telefax +49 (0)781 57211

[info@huber-online.com](mailto:info@huber-online.com)  
[www.huber-online.com](http://www.huber-online.com)

Technischer Service: +49 (0)781 9603-244

-125 °C ... +425 °C

**huber**