



| Elmasonic xtra TT | | | | |
|--------------------------------|------------|------------|-------------|-------------|
| | 30H | 60H | 120H | 200H |
| 220-240 VAC Item No. | 107 6475 | 107 6476 | 107 6477 | 107 6478 |
| 115-120 VAC Item No. | 107 6479 | 107 6480 | 107 6481 | 107 6482 |

Description of Elmasonic xtra TT

Elmasonic xtra TT tabletop ultrasonic devices are designed for use in production environments, workshops and servicing. Features like the permanently integrated sweep function, the switchable dynamic function or the individually settable limit temperature with LED warning indicator make the cleaning of parts in day to day work much easier. All functions and cleaning parameters like time and heating preselection and ultrasonic frequency are arranged user-friendly on the control panel on the front of the device where users can set and monitor them at a glance.




With their ultrasonic performance of 37 kHz, Elmasonic xtra TT devices are rated for continuous heavy-duty operation. The tanks are made of special stainless steel, which makes them robust with a long service life. If they are used in single-shift operation, they are guaranteed for three years.



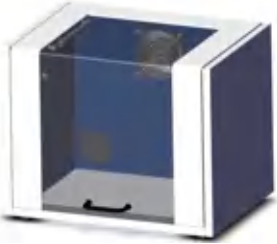
Product features of Elmasonic xtra TT

| | 30H | 60H | 120H | 200H |
|---|------------|------------|-------------|-------------|
| • ultrasonic cleaning device with stainless steel cover | • | • | • | • |
| • efficient 37 kHz ultrasonic high-performance transducers | • | • | • | • |
| • cleaning tank made of cavitation-resistant stainless steel | • | • | • | • |
| • user-friendly and clearly arranged operating panel; sealed against liquid intrusion to protect the electronics | • | • | • | • |
| • LED display for ultrasonic function | • | • | • | • |
| • rotary switch for easy preselection of the cleaning line | • | • | • | • |
| • continuous operating or short-term operation between 1 and 30 minutes | • | • | • | • |
| • permanent sweep function for optimized sound field distribution through frequency modulation | • | • | • | • |
| • Switchable dynamic function increases the peak ultrasonic performance. This improves the effectivity of cleaning and makes it possible to remove even the most stubborn contaminants. | • | • | • | • |
| • dry-run protected heating | • | • | • | • |
| • temperature selectable by rotary switch; from 30°C to 80°C in increments of 5°C (for H units) | • | • | • | • |
| • LED display for heating operation (for H units) | • | • | • | • |
| • detachable power cord with IEC plug | • | • | • | • |
| • selectable temperature threshold with alarm | • | • | • | • |
| • automatic safety shutdown after 8 hours of continuous operation | • | • | • | • |
| • automatic safety shutdown at a bath temperature of 90°C | • | • | • | • |

| Technical data: | | | | |
|---|---|--------------------|--------------------|--------------------|
| Elmasonic xtra TT | 30H | 60H | 120H | 200H |
| Mechanical data | | | | |
| Max. volume of trough (L) | 3 | 6,5 | 14 | 18 |
| Working volume of trough (L) | 2 | 5 | 10 | 13 |
| Internal dimensions of trough W x D x H (approximate in mm) | 240x130x100 | 300x150x150 | 300x240x200 | 320x280x200 |
| External dimensions of device W x D x H (approximate in mm) | 360x230x250 | 420x250x300 | 440x340x350 | 420x390x330 |
| Weight (kg) | 6 | 10 | 11 | 13 |
| Max. basket loading (kg) | 1 | 3 | 7 | 8 |
| Ball valve (") | 3/8 | 3/8 | 3/8 | 3/8 |
| Electrical data | | | | |
| Mains voltage (VAC) | 115-120 220-240 | 115-120 220-240 | 115-120 220-240 | 115-120 220-240 |
| Ultrasound frequency (kHz) | 37 | 37 | 37 | 37 |
| Total power consumption (W) | 540 | 550 | 1000 | 1500 |
| Effective ultrasound output (W) | 140 | 150 | 200 | 300 |
| Max. peak ultrasound output * (W) | 560 | 600 | 800 | 1200 |
| Heating power (W) | 400 | 400 | 800 | 1200 |
| * Due to the signal shape, four times the value results for the maximum peak value of ultrasound output | | | | |
| Noise level | | | | |
| Sound pressure level (L _{pAU})** | <80 dB ** Measured sound pressure level with basket and cover at a distance of 1 metre | | | |
| Ultrasound level (L _{pZ})** | <110 dB ** Measured sound pressure level with basket and cover at a distance of 1 metre | | | |

Accessories:

| Elmasonic xtra TT | | 30H | 60H | 120H | 200H |
|--|---|--|--|---|---|
|  | Stainless steel basket (W x D x H mm): Mesh (mm): Item No. | 196x112x39 6x1 108 2417 | 257x132x74 6x1 108 2408 | 262x222x114 6x1 108 2422 | 280x250x114 6x1 108 2392 |
| | | | | | |
|  | Cover inclusive rubber ring ø 90mm or ø 76mm (with ring) Item No. | 100 4535 | 100 4536 | --- | --- |
| | | | | | |
|  | Glass beaker with lid 600ml, ø 90mm Item No. | 104 6009 | --- | --- | --- |
| | Glass beaker with lid 1000ml, ø 90mm Item No. | 104 6010 | --- | --- | --- |

| | | | | |
|--|---|-----------------------------------|-----|-----|
|  | Glass beaker without lid 600ml, ø 90mm Item No. | 207 034 1000 | --- | --- |
| | Glass beaker without lid 1000ml, ø 90mm Item No. | 207 035 1000 | --- | --- |
|  | Separate lid for glass beakers (ø 90mm, PP blue) Item No. | 104 5677 | --- | --- |
|  | Noise protection box for Elmasonic xtra TT 30 H Item No. | 220 – 240V 200 000 0368 | --- | --- |
| | Item No. | 115 – 120V 103 8838 | --- | --- |

We reserve the right to make technical and visual modifications 11/2018